



Update _2_0_225

Clinical Reminders

**VA-OSP TH INTAKE ASSESMENT UPDATE
Install Guide**

August 2021

Product Development
Office of Information Technology
Department of Veterans Affairs

Table of Contents

INTRODUCTION 1
INSTALL DETAILS 4
INSTALL EXAMPLE 4
POST INSTALLATION 7

Introduction

The original Suicide Prevention Telehealth Interfacility Consult (IFC) and Intake Assessment template were developed by the Veterans Health Administration Office of Mental Health and Suicide Prevention (OMHSP) in collaboration with subject matter experts within the SP 2.0 Clinical Telehealth Workgroup and nationally released earlier in FY21. These tools were designed to create a referral system within the Veterans Integrated Service Network (VISN) Clinical Resource Hubs (CRH) for Veterans with a history of suicidal self-directed violence to receive suicide-prevention psychotherapy interventions via telehealth. A focus group of therapists across several VISN CRHs, who have been using the original intake assessment template, provided feedback about improvements to the tool. Additionally, the SP 2.0 Clinical Telehealth program has learned lessons from both the Human Factors Engineering (HFE) workgroup review and the pilot process for Cognitive Behavior Therapy (CBT-SP) note templates, and these lessons learned have informed our edits to improve the quality of the intake assessment template and have contributed to our next version of this tool.

Referring providers (LIPs) located in Spoke and Hub telehealth facilities from general mental health, specialty mental health, residential mental health and inpatient mental health will complete the referring IFC. Suicide Prevention Telehealth psychiatrists, psychologists, social workers and licensed professional mental health counselors located in the Hub site of the Clinical Resource Hubs will receive and respond to the consults using the Intake Assessment Template. The purpose of the templates is to use health factors to capture data, and collect information about the referral sources and triage responses.

UPDATE_2_0_225 contains 1 Reminder Exchange entry:

UPDATE_2_0_225 VA-OSP TH INTAKE ASSESMENT UPDATE

The exchange file contains the following components:

REMINDER GENERAL FINDINGS

VIEW PROGRESS NOTE TEXT

TIU TEMPLATE FIELD

WP 2/50 NOT REQ

WP 3/70

VA-OSP SP RISK MGMT

VA-OSP SRM CONSULT EMAIL

DATE

EB 20/240

DISPLAY TEXT SP OTHER LOCATION

EB 20/240 REQ

VA-HH:MM

DATE (REQ)

HEALTH FACTORS

VA-REMINDER UPDATES [C]
VA-UPDATE_2_0_225
VA-SUICIDE SAFETY PLAN [C]
VA-SP SAFETY PLAN REVIEWED
VA-OSP TH INTAKE ASSESSMENT [C]
VA-OSP TH SAFETY PLAN UPDATED
VA-OSP TH SAFETY PLAN NEW
VA-OSP TH NOT REC OTHER
VA-OSP TH NOT REC NOT INDICATED
VA-OSP TH NOT REC INCLUSION
VA-OSP TH NOT REC RECEIVING TX
VA-OSP TH NOT REC DIFFERENT ISSUE
VA-OSP TH PST-SP CONSENTED
VA-OSP TH DBT CONSENTED
VA-OSP TH CBT-SP CONSENTED
VA-OSP TH SAFETY PLANNING CONSENTED
VA-OSP TH DECLINED OTHER
VA-OSP TH DECLINED NOT MAIN CONCERN
VA-OSP TH DECLINED DISCUSS SUICIDE
VA-OSP TH DECLINED VIRTUAL
VA-OSP TH DECLINED MOTIVATION
VA-OSP TH DECLINED ATTENDANCE
VA-OSP TH TREATMENT DISCUSSION NONE
VA-OSP TH PST-SP DISCUSSED
VA-OSP TH DBT DISCUSSED
VA-OSP TH CBT-SP DISCUSSED
VA-OSP TH SAFETY PLANNING DISCUSSED
VA-SUICIDE RISK MANAGEMENT [C]
VA-SRM FU LOW ACUTE RISK
VA-SRM FU INTERMEDIATE ACUTE RISK
VA-SRM FU HIGH ACUTE RISK
VA-SUICIDE RISK ASSESSMENT [C]
VA-SUICIDE RISK ACUTE LOW
VA-SUICIDE RISK ACUTE INTERMEDIATE
VA-SUICIDE RISK ACUTE HIGH
VA-OSP TH DISCHARGE [C]
VA-OSP TH ACUTE RISK LOW
VA-OSP TH ACUTE RISK INTERMEDIATE
VA-OSP TH ACUTE RISK HIGH
VA-OSP TH CHRONIC RISK LOW
VA-OSP TH CHRONIC RISK INTERMEDIATE
VA-OSP TH CHRONIC RISK HIGH
VA-SUICIDE RISK EVALUATION TYPE [C]
VA-SUICIDE RISK UPDATE EVALUATION
VA-SUICIDE RISK NEW EVALUATION

VA-OSP TH FACE TO FACE
VA-OSP TH TELEPHONE
VA-OSP TH SESSION PARTIALLY VIDEO
VA-OSP TH VIDEO TO ATLAS
VA-OSP TH VIDEO TO CLINIC
VA-OSP TH VA DEVICE NO
VA-OSP TH VA DEVICE YES
VA-OSP TH VIDEO TO HOME OUTSIDE
VA-OSP TH VIDEO TO HOME RESIDENCE
VA-OSP TH REFERRAL OTHER
VA-OSP TH REFERRAL MH SPECIALTY
VA-OSP TH REFERRAL INPATIENT
VA-OSP TH REFERRAL RRTP
VA-OSP TH REFERRAL SP TEAM
VA-OSP TH REFERRAL PCMHI
VA-OSP TH REFERRAL MHC/BHIP
VA-OSP TH INTAKE ASSESSMENT NOTE

REMINDER SPONSOR

VHA OFFICE OF MENTAL HEALTH AND SUICIDE PREVENTION

REMINDER TERM

VA-REMINDER UPDATE_2_0_225

HEALTH SUMMARY TYPE

VA-OSP TH ACUTE RISK LEVELS 6M
VA-CSRE MOST RECENT

HEALTH SUMMARY OBJECTS

VA-OSP TH ACUTE RISK LEVELS 6M
VA-CSRE MOST RECENT

TIU DOCUMENT DEFINITION

VA-OSP TH ACUTE RISK LEVELS 6M
VA-CSRE MOST RECENT

REMINDER DIALOG

VA-OSP TH INTAKE

Install Details

This update is being distributed as a web host file. The address for the host file is:

https://vaww.va.gov/reminders/docs/UPDATE_2_0_225.PRD


The file will be installed using Reminder Exchange, programmer access is not required.

Installation:

This update can be loaded with users on the system. Installation will take less than **5** minutes.

Install Example

To Load the Web Host File. Navigate to Reminder exchange in Vista

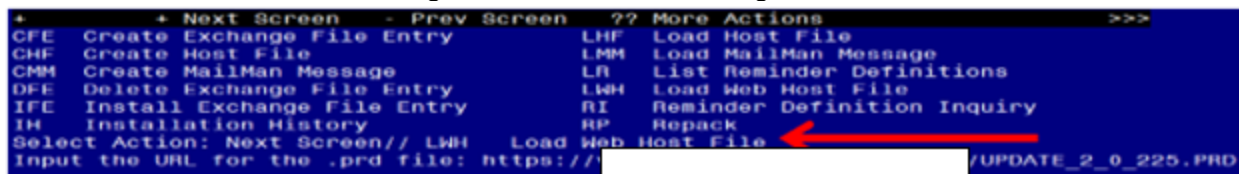

```
+ Next Screen - Prev Screen ?? More Actions >>>
CFE Create Exchange File Entry LHF Load Host File
CHF Create Host File LMM Load MailMan Message
CMM Create MailMan Message LR List Reminder Definitions
DFE Delete Exchange File Entry LWH Load Web Host File
IFE Install Exchange File Entry RI Reminder Definition Inquiry
IH Installation History RP Repack
Select Action: Next Screen// LWH Load Web Host File
Input the URL for the .prd file: https://[redacted]/UPDATE_2_0_225.PRD
```

At the **Select Action:** prompt, enter **LWH** for Load Web Host File


At the **Input the url for the .prd file:** prompt, type the following web address:

https://vaww.va.gov/reminders/docs/UPDATE_2_0_225.PRD

You should see a message at the top of your screen that the file successfully loaded.

[https://\[redacted\]/UPDATE_2_0_225.PRD](https://[redacted]/UPDATE_2_0_225.PRD) successfully loaded.

Search and locate an entry titled UPDATE_2_0_225 VA-OSP TH INTAKE ASSESMENT UPDATE in reminder exchange.

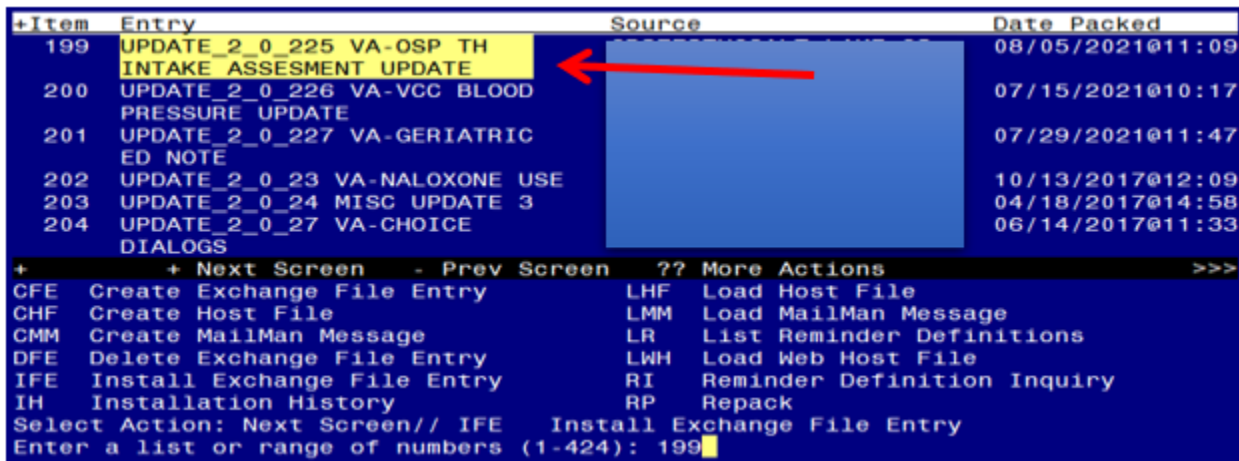

+Item	Entry	Source	Date Packed
199	UPDATE_2_0_225 VA-OSP TH INTAKE ASSESMENT UPDATE	[redacted]	08/05/2021@11:09
200	UPDATE_2_0_226 VA-VCC BLOOD PRESSURE UPDATE	[redacted]	07/15/2021@10:17
201	UPDATE_2_0_227 VA-GERIATRIC ED NOTE	[redacted]	07/29/2021@11:47
202	UPDATE_2_0_23 VA-NALOXONE USE	[redacted]	10/13/2017@12:09
203	UPDATE_2_0_24 MISC UPDATE 3	[redacted]	04/18/2017@14:58
204	UPDATE_2_0_27 VA-CHOICE DIALOGS	[redacted]	06/14/2017@11:33

```
+ Next Screen - Prev Screen ?? More Actions >>>
CFE Create Exchange File Entry LHF Load Host File
CHF Create Host File LMM Load MailMan Message
CMM Create MailMan Message LR List Reminder Definitions
DFE Delete Exchange File Entry LWH Load Web Host File
IFE Install Exchange File Entry RI Reminder Definition Inquiry
IH Installation History RP Repack
Select Action: Next Screen// IFE Install Exchange File Entry
Enter a list or range of numbers (1-424): 199
```

At the **Select Action** prompt, enter **IFE** for Install Exchange File Entry

Enter the number that corresponds with your entry UPDATE_2_0_225 VA-OSP TH INTAKE ASSESMENT UPDATE (in this example it is entry 199. It will vary by site). The date of the exchange file should be 08/05/2021.

```

Component                                     Category  Exists
Source:
Date Packed: 08/05/2021@11:09:44
Package Version: 2.0P42

Description:
The following Clinical Reminder items were selected for packing:
REMINDER DIALOG
  VA-OSP TH INTAKE

REMINDER TERM
  VA-REMINDER UPDATE_2_0_225
Non-exchangeable TIU object(s):
-----
TIU Object: VA-CSRE MOST RECENT NOTE
Object Method: S X=$$TIU^GMTSOBJ(DFN,6600590)

+ Enter ?? for more actions >>>
IA  Install all Components                IS  Install Selected Component
Select Action: Next Screen// IA

```

At the **Select Action** prompt, type **IA** for Install all Components and hit enter.
 Select Action: Next Screen// **IA Install all Components**
 You will see several prompts, for all new entries you will choose **I to Install**
 For components that already exists but the packed component is different, choose **S to Skip**.

You will be prompted to install the reminder dialog component (VA-OSP TH INTAKE). You will choose to install all even though some components may have an X to indicate that they exist.

```

Packed reminder dialog: VA-OSP TH INTAKE [NATIONAL DIALOG]

```

Item	Seq.	Dialog Findings	Type	Exists
1		VA-OSP TH INTAKE	dialog	X
2	10	VA-GP OSP TH INTAKE ASSESSMENT	group	X
		Finding: VA-OSP TH INTAKE ASSESSMENT NOTE (HEALTH FACTOR)		X
3	10.5	VA-OSP TH INTAKE ASSESSMENT VERSION	element	
		Finding: *NONE*		
4	10.10	VA-OSP TH INTAKE DATE TIME	element	X
		Finding: *NONE*		
5	10.20	VA-GP OSP TH INTAKE CHART REVIEW OUTLINE	group	
		Finding: *NONE*		
6	10.20.10	VA-GP OSP TH INTAKE CHART REVIEW	group	X
		Finding: *NONE*		
7	10.20.10.5	VA-BLANK LINE ALT PN TEXT ONLY	element	X
		Finding: *NONE*		
8	10.30	VA-GP OSP TH INTAKE VERIFY ID OUTLINE	group	

```

+ + Next Screen - Prev Screen ?? More Actions
DD Dialog Details DT Dialog Text IS Install Selected
DF Dialog Findings DU Dialog Usage QU Quit
DS Dialog Summary IA Install All
Select Action: Next Screen// IA

```

At the **Select Action** prompt, type **IA** to install the dialog: **VA-OSP TH INTAKE**
 Select Action: Next Screen// **IA Install All**
 Install reminder dialog and all components with no further changes: Y// **Yes**

```

Packed reminder dialog: VA-OSP TH INTAKE [NATIONAL DIALOG]
VA-OSP TH INTAKE (reminder dialog) installed from exchange file.
Item  Seq. Dialog Findings                                     Type  Exists
  1    VA-OSP TH INTAKE                                       dialog  X
  2    10 VA-GP OSP TH INTAKE ASSESSMENT                       group  X
      Finding: VA-OSP TH INTAKE ASSESSMENT NOTE (HEALTH FACTOR)  X
  3    10.5 VA-OSP TH INTAKE ASSESSMENT VERSION               element  X
      Finding: *NONE*
  4    10.10 VA-OSP TH INTAKE DATE TIME                       element  X
      Finding: *NONE*
  5    10.20 VA-GP OSP TH INTAKE CHART REVIEW OUTLINE         group  X
      Finding: *NONE*
  6    10.20.10 VA-GP OSP TH INTAKE CHART REVIEW              group  X
      Finding: *NONE*
  7    10.20.10.5 VA-BLANK LINE ALT PN TEXT ONLY              element  X
      Finding: *NONE*
  8    10.30 VA-GP OSP TH INTAKE VERIFY ID OUTLINE            group  X
+      + Next Screen  - Prev Screen  ?? More Actions
DD  Dialog Details      DT  Dialog Text      IS  Install Selected
DF  Dialog Findings    DU  Dialog Usage     QU  Quit
DS  Dialog Summary     IA  Install All
Select Action: Next Screen// Q

```

After completing this dialog install, you will type **Q**.

```

Component                                     Category  Exists
Source:
Date Packed: 08/05/2021@11:09:44
Package Version: 2.0P42

Description:
The following Clinical Reminder items were selected for packing:
REMINDER DIALOG
  VA-OSP TH INTAKE

REMINDER TERM
  VA-REMINDER UPDATE_2_0_225
Non-exchangeable TIU object(s):
-----
TIU Object: VA-CSRE MOST RECENT NOTE
Object Method: S X=$$TIU^GMTSOBJ(DFN,6600590)

+      + Next Screen  - Prev Screen  ?? More Actions  >>>
IA  Install all Components      IS  Install Selected Component
Select Action: Next Screen// Q

```

You will then be returned to this screen. At the **Select Action** prompt, type **Q**.

Install complete.

Post Installation

1. You should already have a note titled, **SUICIDE PREVENTION TELEHEALTH INTAKE ASSESMENT**. Open and confirm that the version now displays as **version 2.2**
2. **Add your local data objects.**
Add your local data objects for patient address and patient phone numbers to the dialog text of the following reminder dialog elements:
 - a. VAL-OSP TH INTAKE EMERGENCY INFO VTH PHONE
 - b. VAL-OSP TH INTAKE EMERGENCY INFO VTH PT LOCATION HOME

Using the menu Reminder Dialog Management, Reminder Dialogs, change your view to Elements.

Dialog List Mar 03, 2021@08:30:32 Page: 1 of 35
REMINDER VIEW (ALL REMINDERS BY NAME)

Item Reminder Name Linked Dialog Name & Dialog Status

+ Enter ?? for more actions >>>
AR All reminders LR Linked Reminders QU Quit
CV Change View RN Name/Print Name
Select Item: Next Screen// **CV** Change View **<Enter>**

Select one of the following:

D Reminder Dialogs
E Dialog Elements
F Forced Values
G Dialog Groups
P Additional Prompts
R Reminders
RG Result Group (Mental Health)
RE Result Element (Mental Health)

TYPE OF VIEW: R// **E <Enter>**

Next, type *SL* for search, press **<Enter>, then search for **VAL-OSP TH INTAKE EMERGENCY INFO**.**

+ + Next Screen - Prev Screen ?? More Actions >>>
AD Add PT List/Print All QU Quit
CO Copy Dialog INQ Inquiry/Print
CV Change View TE Dialog Taxonomy Edit

Select Item: Next Screen// **SL <Enter>**
Search for: **VAL-OSP TH INTAKE EMERGENCY INFO <Enter>**
...searching for 'VAL-OSP TH INTAKE EMERGENCY INFO'

Dialog List Mar 03, 2021@08:49:59 Page: 1973 of 2012
DIALOG VIEW (DIALOG ELEMENTS)

+Item	Dialog Name	Dialog type	Status
27615	VAL-OSP TH INTAKE EMERGENCY INFO VTH P	Dialog Element	
27616	VAL-OSP TH INTAKE EMERGENCY INFO VTH P	Dialog Element	
+	+ Next Screen	- Prev Screen	?? More Actions >>>

Find Next 'VAL-OSP TH INTAKE EMERGENCY INFO'? Yes// **NO <Enter>**

After searching, you will see the sequence number for the elements in the left column. In this example, I am selecting 27615.

Select Item: Next Screen// **27615 <Enter>**

Dialog Name: VAL-OSP TH INTAKE EMERGENCY INFO VTH PHONE

Current dialog element/group name: VAL-OSP TH INTAKE EMERGENCY INFO VTH PHONE

Used by: VA-GP OSP TH INTAKE EMERGENCY INFO VTH PHONE (Dialog Group)

NAME: VAL-OSP TH INTAKE EMERGENCY INFO VTH PHONE

Replace **<Enter>**

DISABLE: **<Enter>**

CLASS: LOCAL// **<Enter>**

SPONSOR: VHA Office of Mental Health and Suicide Prevention

// **<Enter>**

REVIEW DATE: **<Enter>**

RESOLUTION TYPE: **<Enter>**

ORDERABLE ITEM: **<Enter>**

FINDING ITEM: **<Enter>**

Additional findings: none **<Enter>**

Select ADDITIONAL FINDING: **<Enter>**

DIALOG/PROGRESS NOTE TEXT:

Phone: ****INSERT PATIENT PHONE OBJECT****

Edit? NO// **YES <Enter>**

Add your local data object, then repeat the steps as needed to add your local data object for patient address to the remaining element VAL-OSP TH INTAKE EMERGENCY INFO VTH PT LOCATION HOME.