



Update _2_0_181

Clinical Reminders

**VA-COVID-19 POC ANTIGEN RESULTS NOTE
Install Guide**

December 2020

Product Development
Office of Information Technology
Department of Veterans Affairs

Table of Contents

INTRODUCTION1

INSTALL DETAILS2

INSTALL EXAMPLE.....2


POST INSTALLATION.....5

Introduction

DESCRIPTION:

This template is for use by a Licensed Independent Practitioner (LIP) or Advanced Practice Provider (APP) interpreting the COVID-19 POC Antigen test and directing next steps for testing and/or patient instructions.

The only option that will set the CPRS COVID-19 banner is SYMPTOMATIC POSITIVE, it will set the banner to clinically positive if there are no other positive PCR tests in the lab package.

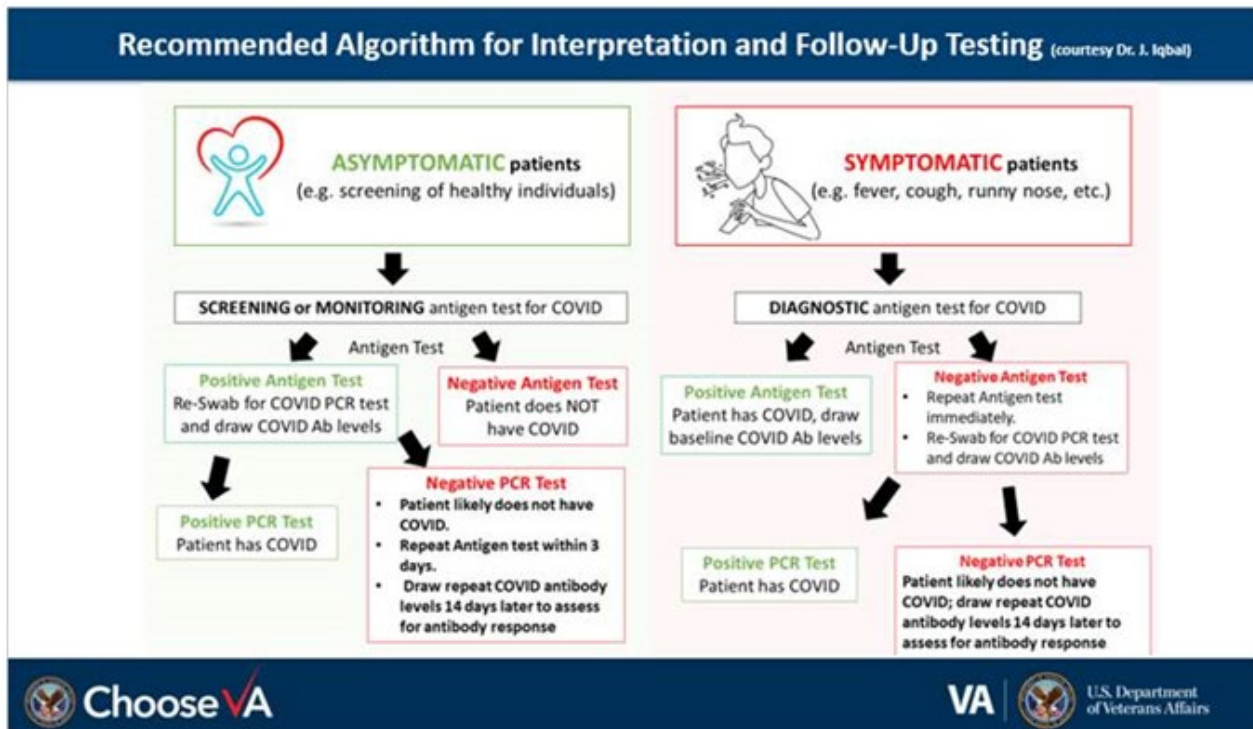

UPDATE_2_0_181 contains 1 Reminder Exchange entry:

UPDATE_2_0_181 VA-COVID-19 POC ANTIGEN RESULTS NOTE

The exchange file contains the following components:

HEALTH FACTORS

VA-REMINDER UPDATES

VA-UPDATE_2_0_181

VA-COVID-19

VA-COVID-19 SUSPECTED

REMINDER TERM

VA-REMINDER UPDATE_2_0_181

REMINDER DIALOG

VA-COVID-19 POC ANTIGEN RESULTS NOTE

Install Details

This update is being distributed as a web host file. The address for the host file is:
https://REDACTED /UPDATE_2_0_181.PRD

The file will be installed using Reminder Exchange, programmer access is not required.

This update can be loaded with users on the system. Installation will take less than 5 minutes.

Install Example

To Load the Web Host File. Navigate to Reminder exchange in Vista

```
+          + Next Screen  - Prev Screen  ?? More Actions  >>>
CFE Create Exchange File Entry          LHF Load Host File
CHF Create Host File                    LMM Load MailMan Message
CMM Create MailMan Message              LR  List Reminder Definitions
DFE Delete Exchange File Entry          LWH Load Web Host File
IFE Install Exchange File Entry         RI  Reminder Definition Inquiry
IH  Installation History                RP  Repack
Select Action: Next Screen// LWH  Load Web Host File
Input the URL for the .prd file: https://[REDACTED] UPDATE_2_0_181.PRD
```

At the **Select Action:** prompt, enter **LWH** for Load Web Host File

At the **Input the url for the .prd file:** prompt, type the following web address:

https://REDACTED /UPDATE_2_0_181.PRD

You should see a message at the top of your screen that the file loaded successfully

```
https://[REDACTED] UPDATE_2_0_181.PRD successfully loaded.
```

Search and locate an entry titled **UPDATE_2_0_181 VA-COVID-19 POC ANTIGEN RESULTS**

+Item	Entry	Source	Date Packed
139	UPDATE_2_0_181 VA-COVID-19 POC ANTIGEN RESULTS		12/22/2020@08:54
140	UPDATE_2_0_182 VA-OSP SUICIDE SAFETY PLAN UPDATE		12/16/2020@12:41
141	UPDATE_2_0_19 VA-CSC INITIAL CAREGIVER ASSESSMENT		01/31/2017@12:14
142	UPDATE_2_0_2 CAREGIVER ELIG		02/18/2016@13:10
143	UPDATE_2_0_20 BREAST CANCER DIALOGS		02/06/2017@12:19
144	UPDATE_2_0_20 COLON CANCER		02/07/2017@09:19
+ + Next Screen - Prev Screen ?? More Actions >>>			
CFE	Create Exchange File Entry	LHF	Load Host File
CHF	Create Host File	LMM	Load MailMan Message
CMM	Create MailMan Message	LR	List Reminder Definitions
DFE	Delete Exchange File Entry	LWH	Load Web Host File
IFE	Install Exchange File Entry	RI	Reminder Definition Inquiry
IH	Installation History	RP	Repack
Select Action: Next Screen// IFE Install Exchange File Entry			
Enter a list or range of numbers (1-368): 139			

At the **Select Action** prompt, enter **IFE** for Install Exchange File Entry
Enter the number that corresponds with your entry **UPDATE_2_0_181 VA-COVID-19 POC ANTIGEN RESULTS**. (in this example it is entry 139 it will vary by site). The date of the exchange file should be 12/22/2020.

Component	Category	Exists
Source: [REDACTED]		
Date Packed: 12/22/2020@08:54:49		
Package Version: 2.0P46		
Description:		
The following Clinical Reminder items were selected for packing:		
REMINDER DIALOG		
VA-COVID-19 POC ANTIGEN RESULTS NOTE		
REMINDER TERM		
VA-REMINDER UPDATE_2_0_181		
Keywords:		
Components:		
+ Enter ?? for more actions >>>		
IA	Install all Components	IS Install Selected Component
Select Action: Next Screen// IA		

At the **Select Action** prompt, type **IA** for Install all Components and hit enter.
Select Action: Next Screen// **IA Install all Components**
You will see several prompts, for all new entries you will choose **I to Install**

You will be prompted to install the reminder dialog component – You will choose to install all even though some components may have an X to indicate that they exist.

```
Packed reminder dialog: VA-COVID-19 POC ANTIGEN RESULTS NOTE [NATIONAL DIALOG]
Item  Seq. Dialog Findings                                     Type  Exists
  1    VA-COVID-19 POC ANTIGEN RESULTS NOTE                   dialog
  2    5 VA-COVID-19 POC ANTIGEN RESULTS NOTE VERSION         element
      Finding: *NONE*
  3    10 VA-GP COVID-19 ANTIGEN TESTING PATIENT TYPE        group
      Finding: *NONE*
  4    10.5 VA-GP COVID-19 ANTIGEN TESTING SYMPTOMATIC        group
      PATIENT
      Finding: *NONE*
  5    10.5.10 VA-GP COVID-19 ANTIGEN TESTING SYMPTOMATIC    group
      PATIENT POS
      Finding: VA-COVID-19 SUSPECTED (HEALTH FACTOR)
  6    10.5.10.5 VA-COVID-19 ANTIGEN TESTING SYMPTOMATIC     element
      PATIENT POS
+      + Next Screen  - Prev Screen  ?? More Actions
DD  Dialog Details      DT  Dialog Text      IS  Install Selected
DF  Dialog Findings    DU  Dialog Usage    QU  Quit
DS  Dialog Summary     IA  Install All
Select Action: Next Screen// IA
```

At the **Select Action** prompt, type **IA** to install the dialog – **VA-COVID-19 POC ANTIGEN RESULTS NOTE**

Select Action: Next Screen// **IA Install All**

Install reminder dialog and all components with no further changes: Y// **Yes**

```
Packed reminder dialog: VA-COVID-19 POC ANTIGEN RESULTS NOTE [NATIONAL DIALOG]
VA-COVID-19 POC ANTIGEN RESULTS NOTE (reminder dialog) installed from exchange f
Item  Seq. Dialog Findings                                     Type  Exists
  1    VA-COVID-19 POC ANTIGEN RESULTS NOTE                   dialog  X
  2    5 VA-COVID-19 POC ANTIGEN RESULTS NOTE VERSION         element  X
      Finding: *NONE*
  3    10 VA-GP COVID-19 ANTIGEN TESTING PATIENT TYPE        group    X
      Finding: *NONE*
  4    10.5 VA-GP COVID-19 ANTIGEN TESTING SYMPTOMATIC        group    X
      PATIENT
      Finding: *NONE*
  5    10.5.10 VA-GP COVID-19 ANTIGEN TESTING SYMPTOMATIC    group    X
      PATIENT POS
      Finding: VA-COVID-19 SUSPECTED (HEALTH FACTOR)
  6    10.5.10.5 VA-COVID-19 ANTIGEN TESTING SYMPTOMATIC     element  X
      PATIENT POS
+      + Next Screen  - Prev Screen  ?? More Actions
DD  Dialog Details      DT  Dialog Text      IS  Install Selected
DF  Dialog Findings    DU  Dialog Usage    QU  Quit
DS  Dialog Summary     IA  Install All
Select Action: Next Screen// Q
```

After completing this dialog install, you will type Q.

Component	Category	Exists
Source: [REDACTED]		
Date Packed: 12/22/2020@08:54:49		
Package Version: 2.0P46		
Description:		
The following Clinical Reminder items were selected for packing:		
REMINDER DIALOG		
VA-COVID-19 POC ANTIGEN RESULTS NOTE		
REMINDER TERM		
VA-REMINDER UPDATE_2_0_181		
Keywords:		
Components:		
+ Enter ?? for more actions >>>		
IA	Install all Components	IS Install Selected Component
Select Action: Next Screen// Q		

You will then be returned to this screen. At the **Select Action** prompt, type **Q**.

Install complete.

Post Installation

1. Setup of Note title

You will need to create a new note title for this template. The note title needed is:
COVID-19 POC ANTIGEN RESULTS NOTE

The note title should be mapped to the following VHA Enterprise Standard Title:
 INFECTIOUS DISEASE NOTE

Note: Comprehensive information on Creating Document Definitions can be found beginning on page 47 of the [TIU/ASU Implementation Guide](#)

2. Associate the new reminder dialogs with the note title in CPRS

To do this, in VistA go to Reminder Manager, then CPRS Reminder Configuration then LINK- Link Reminder Dialog to Template.

Type in the following information:

SELECT CPRS REMINDER CONFIGURATION <TEST ACCOUNT> OPTION: **LINK** LINK
 REMINDER DIALOG TO TEMPLATE

SELECT DIALOG DEFINITION: **VA-COVID-19 POC ANTIGEN RESULTS NOTE**
 REMINDER DIALOG NATIONAL
 ...OK? YES// (YES)

ENTER TEMPLATE NAME: COVID-19 POC ANTIGEN RESULTS NOTE

LINK TEMPLATE TO DOCUMENT TITLE? YES

SELECT DOCUMENT DEFINITION: COVID-19 POC ANTIGEN RESULTS NOTE

3. Edit the reminder dialog element VAL-COVID-19 CONTACT VA NURSE HOTLINE and replace “(XXX) XXX-XXXX” in the dialog text with your facility VA nurse hotline phone information.
4. Add your local quick order for COVID-19 PCR testing as a finding item to the reminder dialog element **VAL-COVID-19 ANTIGEN TESTING ORDER**.
5. Enable the COVID-19 Antibody testing order element if appropriate for your site and add your local quick order for COVID-19 Antibody testing as a finding item to the reminder dialog element **VAL-COVID-19 ANTIBODY TESTING ORDER**.

Before enabling the Antibody test in this template, coordinate with your Lab Manager and Infection Control/Infectious Disease staff.

If it is enabled, consider editing the reminder dialog text and adding directions for various testing locations to ensure COVID-19 positive patients do not have to move about the facility to get tested. E.g. AntiBODY testing can be added in the ED and the COVID-19 Pod where blood draw is available.

The order is used in multiple places in the dialog, when enabled, it will show as an ordering option:

Example 1:

POC COVID-19 Antigen testing done for:

Symptomatic patient

POSITIVE result

++Patient notified this is considered definitive testing.++

COVID-19 Antibody testing ordered

Example 2:

Screening (Pre-appointment/Pre-procedure testing)

POSITIVE result

++Patient notified this is NOT considered definitive testing.++

COVID-19 PCR testing ordered to confirm

COVID-19 Antibody testing ordered

6. Add local content if needed - OPTIONAL

This template was built with the standard follow up instructions included. If a site would like to add guidance or information specific to that site, you can add text to the local reminder dialog element VAL-COVID-19 ANTIGEN TESTING POS LOCAL INFO for patients with a positive result and VAL-COVID-19 ANTIGEN TESTING SYMPTOMATIC PATIENT NEG for symptomatic patients with negative POC results. These dialog elements are at the end of the sections. This dialog elements will install disabled and you will need to remove the disable flag for them to be used.