



**Update \_2\_0\_223**

**Clinical Reminders**

**VA-REMOTE TEMPERATURE MONITORING  
Install Guide**

**September 2021**

Product Development  
Office of Information Technology  
Department of Veterans Affairs

**Table of Contents**

**INTRODUCTION ..... 1**  
**INSTALL DETAILS ..... 1**  
**INSTALL EXAMPLE ..... 2**  
**POST INSTALLATION ..... 5**

## Introduction

---

The Remote Temperature Monitoring Template should be used at any site with patients enrolled in Remote Temperature Monitoring. Separate enrollment notes should be entered for any patient that is enrolled in the program. Escalation alert notes should be entered when escalation alerts are received from the vendor via email. Disenrollment notes should be entered when a patient is disenrolled from the program for any reason. These notes will include health factors which will make it much easier to track who has been issued a smart mat nationally as well as track outcomes for this program.

Any sites using this program should create a REMOTE TEMPERATURE MONITORING CLINIC with the primary stop code of the service (example podiatry, PT, etc) and the secondary stop code of 698 for Remote Monitoring. Each of these notes will be entered into this clinic to capture workload. Coding is explained in the above flowchart.

UPDATE\_2\_0\_223 contains 1 Reminder Exchange entry:

**UPDATE\_2\_0\_223 VA-REMOTE TEMPERATURE MONITORING**

**The exchange file contains the following components:**

### **TIU TEMPLATE FIELD**

TEXT (1-10 CHAR) REQ

WP FREE TEXT (74X4)

TEXT (1-60 CHAR)R

### **HEALTH FACTORS**

VA-REMINDER UPDATES [C]

VA-REMOTE TEMPERATURE MONITORING [C]

VA-UPDATE\_2\_0\_223

VA-REMOTE TEMPERATURE MONITORING-ALERT

VA-REMOTE TEMPERATURE MONITORING-DISENROLL

VA-REMOTE TEMPERATURE MONITORING-ENROLL

### **REMINDER TERM**

VA-REMINDER UPDATE\_2\_0\_223

### **REMINDER DIALOG**

VA-REMOTE TEMPERATURE MONITORING

### **REMINDER SPONSOR**

OFFICE OF REHABILITATION SERVICES IN VACO

# Install Details

---

This update is being distributed as a web host file. The address for the host file is:

**https://REDACTED /UPDATE\_2\_0\_223.PRD**

The file will be installed using Reminder Exchange, programmer access is not required.

Installation:

=====

This update can be loaded with users on the system. Installation will take less than **5** minutes.

# Install Example

---

To Load the Web Host File. Navigate to Reminder exchange in Vista

```
CFE Create Exchange File Entry      LHF Load Host File
CHF Create Host File                LMM Load MailMan Message
CMM Create MailMan Message          LR List Reminder Definitions
DFE Delete Exchange File Entry      LWH Load Web Host File
IFE Install Exchange File Entry      RI Reminder Definition Inquiry
IH Installation History              RP Repack
Select Action: Next Screen// LWH Load Web Host File
Input the URL for the .prd file: https://[REDACTED]/UPDATE_2_0_223_1.PRD
Added Reminder Exchange entry UPDATE_2_0_223 VA-REMOTE TEMPERATURE MONITORING
```

At the **Select Action:** prompt, enter **LWH** for Load Web Host File

At the **Input the url for the .prd file:** prompt, type the following web address:

[https://Redacted/UPDATE\\_2\\_0\\_223.PRD](https://Redacted/UPDATE_2_0_223.PRD)

You should see a message at the top of your screen that the file successfully loaded.

Search and locate an entry titled **UPDATE\_2\_0\_223 VA-REMOTE TEMPERATURE MONITORING** in reminder exchange.

```
+Item Entry Source Date Packed
300 UPDATE_2_0_223 VA-REMOTE TEMPERATURE MONITORING 07/27/2021@14:32
301 UPDATE_2_0_224 VA-DBT AND DISCHARGE SUMMARY 07/15/2021@11:26
302 UPDATE_2_0_224 VA-DBT AND DISCHARGE SUMMARY 07/27/2021@12:44
303 UPDATE_2_0_23 VA-NALOXONE USE 10/13/2017@12:09
304 UPDATE_2_0_24 MISC UPDATE 3 05/08/2017@04:44
305 UPDATE_2_0_26 VA-HEAD AND NECK CANCER DIALOGS 08/22/2017@13:17
+ + Next Screen - Prev Screen ?? More Actions >>>
CFE Create Exchange File Entry LHF Load Host File
CHF Create Host File LMM Load MailMan Message
CMM Create MailMan Message LR List Reminder Definitions
DFE Delete Exchange File Entry LWH Load Web Host File
IFE Install Exchange File Entry RI Reminder Definition Inquiry
IH Installation History RP Repack
Select Action: Next Screen// IFE Install Exchange File Entry
Enter a list or range of numbers (1-543): 300
```

At the **Select Action** prompt, enter **IFE** for Install Exchange File Entry

Enter the number that corresponds with your entry UPDATE\_2\_0\_223 VA-REMOTE TEMPERATURE MONITORING (*in this example it is entry 300. It will vary by site*). The date of the exchange file should be 07/27/2021.

```
Date Packed: 07/27/2021@14:32:53
Package Version: 2.0P42

Description:
The following Clinical Reminder items were selected for packing:
REMINDER DIALOG
  VA-REMOTE TEMPERATURE MONITORING

REMINDER TERM
  VA-REMINDER UPDATE_2_0_223

Keywords:

Components:

+      Enter ?? for more actions      >>>
IA  Install all Components              IS  Install Selected Component
Select Action: Next Screen// IA
```

At the **Select Action** prompt, type **IA** for Install all Components and hit enter.

Select Action: Next Screen// **IA Install all Components**

You will see several prompts, for all new entries you will choose **I to Install**

**Install** all new items.

You will be prompted to install the reminder dialog component (**VA-REMOTE TEMPERATURE MONITORING**). You will choose to install all even though some components may have an X to indicate that they exist.

Item	Seq.	Dialog Findings	Type	Exists		
1		VA-REMOTE TEMPERATURE MONITORING	dialog			
2	5	VA-REMOTE TEMPERATURE MONITORING VERSION Finding: *NONE*	element			
3	10	VA-GP REMOTE THERMAL MONITORING CONNECTION Finding: *NONE*	group			
4	10.5	VA-GP REMOTE THERMAL MONITORING CONNECTION-ENROLLMENT Finding: *NONE*	group			
5	10.5.5	VA-TXT REMOTE TEMPERATURE MONITORING ENROLLMENT-INFO BUTTON Finding: VIEW PROGRESS NOTE TEXT (REMINDER GENERAL FINDING)	element	X		
6	PXRM GF	VIEW BUTTON	prompt	X		
7	10.5.10	VA-HF REMOTE TEMPERATURE	element			
+      + Next Screen      - Prev Screen      ?? More Actions						
DD		Dialog Details	DT	Dialog Text	IS	Install Selected
DF		Dialog Findings	DU	Dialog Usage	QU	Quit
DS		Dialog Summary	IA	Install All		
Select Action: Next Screen// IA						

At the **Select Action** prompt, type **IA** to install the dialog: **VA-REMOTE TEMPERATURE MONITORING**

Select Action: Next Screen// **IA Install All**

Install reminder dialog and all components with no further changes: Y// **Yes**

After completing this dialog install, you will type **Q**.

Item	Seq.	Dialog Findings	Type	Exists		
1		VA-REMOTE TEMPERATURE MONITORING	dialog	X		
2	5	VA-REMOTE TEMPERATURE MONITORING VERSION Finding: *NONE*	element	X		
3	10	VA-GP REMOTE THERMAL MONITORING CONNECTION Finding: *NONE*	group	X		
4	10.5	VA-GP REMOTE THERMAL MONITORING CONNECTION-ENROLLMENT Finding: *NONE*	group	X		
5	10.5.5	VA-TXT REMOTE TEMPERATURE MONITORING ENROLLMENT-INFO BUTTON Finding: VIEW PROGRESS NOTE TEXT (REMINDER GENERAL FINDING)	element	X		
6	PXRM GF	VIEW BUTTON	prompt	X		
7	10.5.10	VA-HF REMOTE TEMPERATURE	element	X		
+      + Next Screen      - Prev Screen      ?? More Actions						
DD		Dialog Details	DT	Dialog Text	IS	Install Selected
DF		Dialog Findings	DU	Dialog Usage	QU	Quit
DS		Dialog Summary	IA	Install All		
Select Action: Next Screen// Q						

```
Date Packed: 07/27/2021@14:32:53
Package Version: 2.0P42
```

```
Description:
The following Clinical Reminder items were selected for packing:
REMINDER DIALOG
  VA-REMOTE TEMPERATURE MONITORING
```

```
REMINDER TERM
  VA-REMINDER UPDATE_2_0_223
```

```
Keywords:
```

```
Components:
```

```
+      + Next Screen  - Prev Screen  ?? More Actions      >>>
IA  Install all Components      IS  Install Selected Component
Select Action: Next Screen// Q
```

You will then be returned to this screen. At the **Select Action** prompt, type **Q**.

Install complete.

## Post Installation

---

### 1. Make the Dialog template available to be attached to a Personal or Shared template or to be able to attach the dialog to a progress note title.

Select Reminder Managers Menu <TEST ACCOUNT> Option: **CP CPRS**

#### Reminder Configuration

CA	Add/Edit Reminder Categories
CL	CPRS Lookup Categories
CS	CPRS Cover Sheet Reminder List
MH	Mental Health Dialogs Active
PN	Progress Note Headers
RA	Reminder GUI Resolution Active
<b>TIU</b>	<b>TIU Template Reminder Dialog Parameter</b>
DEVL	Evaluate Coversheet List on Dialog Finish
DL	Default Outside Location
PT	Position Reminder Text at Cursor
LINK	Link Reminder Dialog to Template
TEST	CPRS Coversheet Time Test
NP	New Reminder Parameters
GEC	GEC Status Check Active
WH	WH Print Now Active

Select CPRS Reminder Configuration <TEST ACCOUNT> Option: **TIU TIU**

#### Template Reminder Dialog Parameter

Reminder Dialogs allowed as Templates may be set for the following:

1	User	USR	[choose from NEW PERSON]
3	Service	SRV	[choose from SERVICE/SECTION]
4	Division	DIV	[SALT LAKE CITY]

5 System SYS [NATREM.FO-SLC.MED.VA.GOV]

Enter selection: **5** Choose the appropriate number for your site. Your site may do this by System or other levels. The example below uses SYSTEM level.

Select Display Sequence: **?**

388 VA-COVID-19 TELEPHONE FOLLOW-UP

389 VA-SRA-CONSULT RESPONSE (D)

When you type a question mark above, you will see the list of #'s (with dialogs) that are already taken. Choose a number **NOT** on this list. For this example, 390 is not present so I will use 390.

Select Display Sequence: **390**

Are you adding 390 as a new Display Sequence? Yes// **YES**

Display Sequence: 390// <Enter>

Clinical Reminder Dialog: **VA-REMOTE TEMPERATURE MONITORING** reminder dialog NATIONAL

...OK? Yes// <Enter> (Yes)

## 2. Setup of Note Title.

The note title to create to associate with this template is: **REMOTE TEMPERATURE MONITORING**


VHA Enterprise Standard Title: **TELEHEALTH NOTE**

*Note: Comprehensive information on Creating new Document Classes and Titles can be found beginning on page 54 of the TIU/ASU Implementation Guide at the following link: [http://www.va.gov/vdl/documents/Clinical/CPRS-Text Integration Utility \(TIU\)/tiuim.pdf](http://www.va.gov/vdl/documents/Clinical/CPRS-Text%20Integration%20Utility%20(TIU)/tiuim.pdf)*

## 3. Associate the reminder dialog with the note title in CPRS

In CPRS, open the Template Editor, expand Document Titles, select the appropriate folder (if applicable), then click New Template (top right-hand corner).



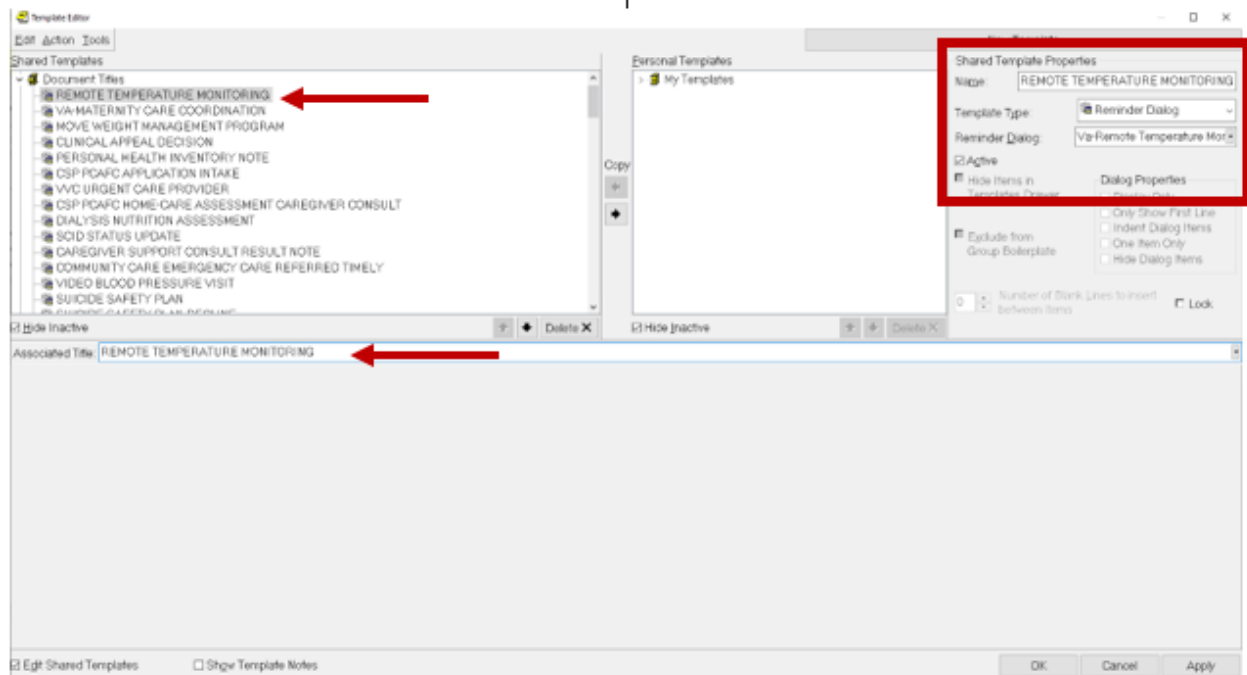

Type in the following information:

Name: **REMOTE TEMPERATURE MONITORING**

Template Type: **Reminder Dialog**

Reminder Dialog: **VA-REMOTE TEMPERATURE MONITORING**

Next, type in the Associated Title name – **REMOTE TEMPERATURE MONITORING** – then click Apply.