



Update_2_0_232

Clinical Reminders

VA-TELEWOUND IMAGER AND READER TEMPLATES

Install Guide

December 2021

Product Development
Office of Information Technology
Department of Veterans Affairs

Table of Contents

INTRODUCTION	1
INSTALL DETAILS	2
INSTALL EXAMPLE	2
POST INSTALLATION	6

Introduction

The imaging templates will be used by Nurses/telehealth care technicians when they image a Veterans wound in the field following a consult placed by Primary care providers/wound care specialist. The settings will vary based on the need. Imaging can take place in outpatient setting, inpatient setting, or into the home. The imaging will also occur in an intra/inter facility capacity. The reader template will be used by wound care specialists following the imaging to assess image quality and provide guidance on wound treatment and management.

UPDATE_2_0_232 contains 1 Reminder Exchange entry:

UPDATE_2_0_232 VA-TELEWOUND IMAGER AND READER TEMPLATES

The exchange file contains the following components:

TIU TEMPLATE FIELD

TEXT (1-50 CHAR) REQ
WP FREE TEXT REQ(74X4)
VA*WP 65/6 INDENT2
TEXT (1-70) REQUIRED
PAIN SCALE 0-10
EDIT 20 (REQ)
EDIT 45(REQUIRED)
WOUND DRESSINGS
EDIT 20
YES NO NA
VA-TELEWOUND BETTER/WORSE/SAME
YES/NO (R/VERT)
YES/NO/NA (R/VERT)

HEALTH FACTORS

VA-REMINDER UPDATES [C]
VA-UPDATE_2_0_232
VA-TELEWOUND IMAGE QUALITY [C]
VA-TWR IMAGE QUALITY NEEDS IMPROVEMENT
VA-TWR IMAGE QUALITY INADEQUATE
VA-TWR IMAGE QUALITY ADEQUATE
VA-TELEWOUND IMAGER TECHNOLOGY [C]
VA-TWI TECHNOLOGY-OTHER
VA-TWI TECHNOLOGY-TE3
VA-TWI TECHNOLOGY-3D IMAGE
VA-TWI TECHNOLOGY-POINT AND SHOOT

REMINDER SPONSOR

OFFICE OF TELEHEALTH SERVICES STORE AND FORWARD TELEHEALTH

REMINDER TERM

VA-REMINDER UPDATE_2_0_232

REMINDER DIALOG

VA-TELEWOUND READER NOTE

VA-TELEWOUND IMAGER NOTE

Install Details

This update is being distributed as a web host file. The address for the host file is:

https://REDACTED/UPDATE_2_0_232.PRD

The file will be installed using Reminder Exchange, programmer access is not required.


Installation:

=====

This update can be loaded with users on the system. Installation will take less than 20 minutes.

Install Example

To Load the Web Host File. Navigate to Reminder exchange in Vista

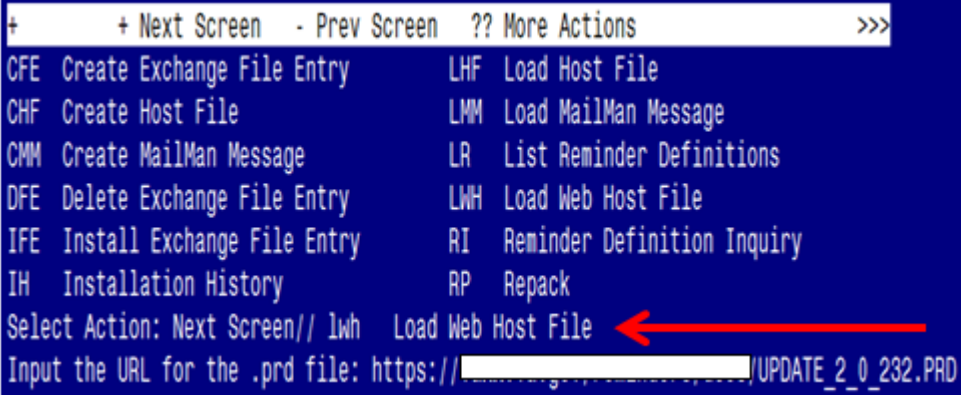

```
+ + Next Screen - Prev Screen ?? More Actions >>>
CFE Create Exchange File Entry LHF Load Host File
CHF Create Host File LMM Load MailMan Message
CMM Create MailMan Message LR List Reminder Definitions
DFE Delete Exchange File Entry LWH Load Web Host File
IFE Install Exchange File Entry RI Reminder Definition Inquiry
IH Installation History RP Repack
Select Action: Next Screen// lwh Load Web Host File
Input the URL for the .prd file: https://REDACTED/UPDATE_2_0_232.PRD
```

At the **Select Action:** prompt, enter **LWH** for Load Web Host File

At the **Input the url for the .prd file:** prompt, type the following web address:

https://REDACTED/UPDATE_2_0_232.PRD

You should see a message at the top of your screen that the file successfully loaded.

https://vaww.va.gov/reminders/docs/UPDATE_2_0_232.PRD successfully loaded.

Search for and locate the reminder exchange entry titled UPDATE_2_0_232 VA-TELEWOUND IMAGER AND READER TEMPLATES.

+Item	Entry	Source	Date Packed
350	UPDATE_2_0_232 VA-TELEWOUND IMAGER AND READER TEMPLATES		12/20/2021@06:37
351	UPDATE_2_0_233 VA-MISC UPDATE 15		10/07/2021@13:11
352	UPDATE_2_0_237 VA-COVID-19 TELEPHONE FOLLOW-UP UPDATE		10/26/2021@04:22
353	UPDATE_2_0_238 VA-PROBLEM SOLVING THERAPY SUICIDE PREVENTION		10/12/2021@07:52
354	UPDATE 2 0 239 VA-CAREGIVER		10/19/2021@12:29

+	Next Screen	- Prev Screen	?? More Actions	>>>
CFE	Create Exchange File Entry	LHF	Load Host File	
CHF	Create Host File	LMM	Load MailMan Message	
CMM	Create MailMan Message	LR	List Reminder Definitions	
DFE	Delete Exchange File Entry	LWH	Load Web Host File	
IFE	Install Exchange File Entry	RI	Reminder Definition Inquiry	
IH	Installation History	RP	Repack	

Select Action: Next Screen// IFE Install Exchange File Entry
Enter a list or range of numbers (1-604): 350

At the **Select Action** prompt, enter **IFE** for Install Exchange File Entry
Enter the number that corresponds with your entry titled **UPDATE_2_0_232 VA-TELEWOUND IMAGER AND READER TEMPLATES** (in this example it is entry 350, it will vary by site). The exchange file date should be 12/20/2021@06:37.

Component	Category	Exists
Source: [REDACTED]		
Date Packed: 12/20/2021@06:37:13		
Package Version: 2.0P71		
Description:		
The following Clinical Reminder items were selected for packing:		
REMINDER DIALOG		
VA-TELEWOUND IMAGER NOTE		
VA-TELEWOUND READER NOTE		
REMINDER TERM		
VA-REMINDER UPDATE_2_0_232		
Keywords:		
Components:		
+ Enter ?? for more actions >>>		
IA	Install all Components	IS Install Selected Component

Select Action: Next Screen// IA

At the **Select Action** prompt, type **IA** for Install all Components and hit enter.
Select Action: Next Screen// IA Install all Components

For any new components choose – **I to INSTALL**

For any existing components choose - **S to SKIP**

You will be prompted to install the 1st dialog component:

```
Packed reminder dialog: VA-TELEWOUND READER NOTE [NATIONAL DIALOG]
+-----+-----+-----+-----+-----+
Item  Seq.  Dialog Findings                                     Type  Exists
-----+-----+-----+-----+-----+
1      1     VA-TELEWOUND READER NOTE                           dialog  X
2      5     VA-TELEWOUND READER NOTE VERSION                     element X
      Finding: *NONE*
3     10     VA-GP TWR WOUND A                                   group
      Finding: *NONE*
4    10.5   VA-TXT TWR WOUND A-LOCATION                           element
      Finding: *NONE*
5    10.10  VA-TXT TWR WOUND A-HISTORY                           element
      Finding: *NONE*
6    10.15  VA-GP TWR WOUND A-RECOMMENDATIONS                     group
      Finding: *NONE*
7    10.15.5 VA-GP TWR WOUND A RECOMMENDATIONS-TREATMENT         group
      Finding: *NONE*
+-----+-----+-----+-----+-----+
+      + Next Screen  - Prev Screen  ?? More Actions
DD  Dialog Details    DT  Dialog Text    IS  Install Selected
DF  Dialog Findings  DU  Dialog Usage   QU  Quit
DS  Dialog Summary   IA  Install All
Select Action: Next Screen// IA
```

At the **Select Action** prompt, type **IA** to install the dialog – **VA-TELEWOUND READER NOTE**

Select Action: Next Screen// **IA Install All**

Install reminder dialog and all components with no further changes: **Y// ES**

```
Packed reminder dialog: VA-TELEWOUND READER NOTE [NATIONAL DIALOG]
VA-TELEWOUND READER NOTE (reminder dialog) installed from exchange file.
+-----+-----+-----+-----+-----+
Item  Seq.  Dialog Findings                                     Type  Exists
-----+-----+-----+-----+-----+
1      1     VA-TELEWOUND READER NOTE                           dialog  X
2      5     VA-TELEWOUND READER NOTE VERSION                     element X
      Finding: *NONE*
3     10     VA-GP TWR WOUND A                                   group   X
      Finding: *NONE*
4    10.5   VA-TXT TWR WOUND A-LOCATION                           element X
      Finding: *NONE*
5    10.10  VA-TXT TWR WOUND A-HISTORY                           element X
      Finding: *NONE*
6    10.15  VA-GP TWR WOUND A-RECOMMENDATIONS                     group   X
      Finding: *NONE*
7    10.15.5 VA-GP TWR WOUND A RECOMMENDATIONS-TREATMENT         group   X
      Finding: *NONE*
+-----+-----+-----+-----+-----+
+      + Next Screen  - Prev Screen  ?? More Actions
DD  Dialog Details    DT  Dialog Text    IS  Install Selected
DF  Dialog Findings  DU  Dialog Usage   QU  Quit
DS  Dialog Summary   IA  Install All
Select Action: Next Screen// Q
```

When the dialog has completed installation, you will then be returned to this screen. At the

Select Action prompt, type **Q**.

You will be prompted to install the 2nd dialog component:

```
Packed reminder dialog: VA-TELEWOUND IMAGER NOTE [NATIONAL DIALOG]
+-----+-----+-----+-----+
Item  Seq. Dialog Findings                               Type  Exists
-----+-----+-----+-----+
  1    VA-TELEWOUND IMAGER NOTE                         dialog  X
-----+-----+-----+-----+
  2    5 VA-TELEWOUND IMAGER NOTE VERSION
       Finding: *NONE*                                 element  X
-----+-----+-----+-----+
  3    10 VA-GP TWI CONSULT REQUEST
       Finding: *NONE*                                 group    X
-----+-----+-----+-----+
  4    10.5 VA-TXT TWI CONSENT/EDUCATION
       Finding: *NONE*                                 element  X
-----+-----+-----+-----+
  5    15 VA-GP TWI TECHNOLOGY USED
       Finding: *NONE*                                 group    X
-----+-----+-----+-----+
  6    15.5 VA-HF TWI TECHNOLOGY USED-POINT/SHOOT
       Finding: VA-TWI TECHNOLOGY-POINT AND SHOOT (HEALTH FACTOR)
                                                                 element  X
-----+-----+-----+-----+
  7    15.10 VA-HF TWI TECHNOLOGY USED-3D IMAGING
                                                                 element  X
-----+-----+-----+-----+
+      + Next Screen  - Prev Screen  ?? More Actions
DD  Dialog Details      DT  Dialog Text      IS  Install Selected
DF  Dialog Findings    DU  Dialog Usage     QU  Quit
DS  Dialog Summary     IA  Install All
Select Action: Next Screen// IA
```

At the **Select Action** prompt, type **IA** to install the dialog – **VA-TELEWOUND IMAGER NOTE**
 Select Action: Next Screen// **IA Install All**
 Install reminder dialog and all components with no further changes: **Y// ES**

```
Packed reminder dialog: VA-TELEWOUND IMAGER NOTE [NATIONAL DIALOG]
VA-TELEWOUND IMAGER NOTE (reminder dialog) installed from exchange file.
+-----+-----+-----+-----+
Item  Seq. Dialog Findings                               Type  Exists
-----+-----+-----+-----+
  1    VA-TELEWOUND IMAGER NOTE                         dialog  X
-----+-----+-----+-----+
  2    5 VA-TELEWOUND IMAGER NOTE VERSION
       Finding: *NONE*                                 element  X
-----+-----+-----+-----+
  3    10 VA-GP TWI CONSULT REQUEST
       Finding: *NONE*                                 group    X
-----+-----+-----+-----+
  4    10.5 VA-TXT TWI CONSENT/EDUCATION
       Finding: *NONE*                                 element  X
-----+-----+-----+-----+
  5    15 VA-GP TWI TECHNOLOGY USED
       Finding: *NONE*                                 group    X
-----+-----+-----+-----+
  6    15.5 VA-HF TWI TECHNOLOGY USED-POINT/SHOOT
       Finding: VA-TWI TECHNOLOGY-POINT AND SHOOT (HEALTH FACTOR)
                                                                 element  X
-----+-----+-----+-----+
  7    15.10 VA-HF TWI TECHNOLOGY USED-3D IMAGING
                                                                 element  X
-----+-----+-----+-----+
+      + Next Screen  - Prev Screen  ?? More Actions
DD  Dialog Details      DT  Dialog Text      IS  Install Selected
DF  Dialog Findings    DU  Dialog Usage     QU  Quit
DS  Dialog Summary     IA  Install All
Select Action: Next Screen// Q
```

When the dialog has completed installation, you will then be returned to this screen. At the **Select Action** prompt, type **Q**.

```

Component                                     Category      Exists
-----
Source:
Date Packed: 12/20/2021@06:37:13
Package Version: 2.0P71

Description:
The following Clinical Reminder items were selected for packing:
REMINDER DIALOG
  VA-TELEWOUND IMAGER NOTE
  VA-TELEWOUND READER NOTE

REMINDER TERM
  VA-REMINDER UPDATE_2_0_232

Keywords:

Components:
+ Next Screen - Prev Screen ?? More Actions >>>
IA  Install all Components          IS  Install Selected Component
Select Action: Next Screen// Q

```

You will be returned to this screen. Type Q to quit.
Install Completed

Post Installation

1. **Make the Dialog templates available to be attached to a Personal or Shared templates or to be able to attach the dialog to a progress note title.**

Select Reminder Managers Menu <TEST ACCOUNT> Option: **CP CPRS Reminder**

Configuration

- CA Add/Edit Reminder Categories
- CL CPRS Lookup Categories
- CS CPRS Cover Sheet Reminder List
- MH Mental Health Dialogs Active
- PN Progress Note Headers
- RA Reminder GUI Resolution Active
- TIU TIU Template Reminder Dialog Parameter**
- DEVL Evaluate Coversheet List on Dialog Finish
- DL Default Outside Location
- PT Position Reminder Text at Cursor
- LINK Link Reminder Dialog to Template
- TEST CPRS Coversheet Time Test
- NP New Reminder Parameters
- GEC GEC Status Check Active
- WH WH Print Now Active

Select CPRS Reminder Configuration <TEST ACCOUNT> Option: **TIU TIU Template Reminder**

Dialog Parameter

Reminder Dialogs allowed as Templates may be set for the following:

- 1 User USR [choose from NEW PERSON]

- 3 Service SRV [choose from SERVICE/SECTION]
- 4 Division DIV [SALT LAKE CITY]
- 5 System SYS [NATREM.FO-SLC.MED.VA.GOV]

Enter selection: **5** Choose the appropriate number for your site. Your site may do this by System or other levels. The example below uses SYSTEM level.

Setting Reminder Dialogs allowed as Templates for System: NATREM.FO-SLC.MED.VA.GOV

Select Display Sequence: **?**

388 VA-COVID-19 TELEPHONE FOLLOW-UP

389 VA-SRA-CONSULT RESPONSE (D)

*When you type a question mark above, you will see the list of #'s (with dialogs) that are already taken. Choose a number **NOT** on this list. For this example, 390 is not present so I will use 390.*

Select Display Sequence: **390**

Are you adding 390 as a new Display Sequence? Yes// **YES**

Display Sequence: 390// **<Enter>** 390

Clinical Reminder Dialog: **VA-TELEWOUND IMAGER NOTE** reminder dialog NATIONAL

...OK? Yes// **<Enter>** (Yes)

REPEAT THIS STEP TO ADD THE 2ND DIALOG AS A SEQUENCE

Display Sequence: 391// **<Enter>** 391

Clinical Reminder Dialog: **VA-TELEWOUND READER NOTE** reminder dialog NATIONAL

...OK? Yes// **<Enter>** (Yes)

2. Setup of Note Titles.

The note titles to create to associate with these templates are:

TELEWOUND IMAGER CONSULT NOTE

TELEWOUND READER CONSULT NOTE


VHA Enterprise Standard Title: **TELEHEALTH CONSULT**

You can add these to your document class where CONSULTS are located.

***Note:** Comprehensive information on Creating new Document Classes and Titles can be found beginning on page 54 of the TIU/ASU Implementation Guide at the following link:*
[http://www.va.gov/vdl/documents/Clinical/CPRS-Text_Integration_UTILITY_\(TIU\)/tiuim.pdf](http://www.va.gov/vdl/documents/Clinical/CPRS-Text_Integration_UTILITY_(TIU)/tiuim.pdf)

3. Associate the reminder dialogs with the note titles in CPRS

In CPRS, open the Template Editor, expand Document Titles, select the appropriate folder (if applicable), then click New Template (top right-hand corner).

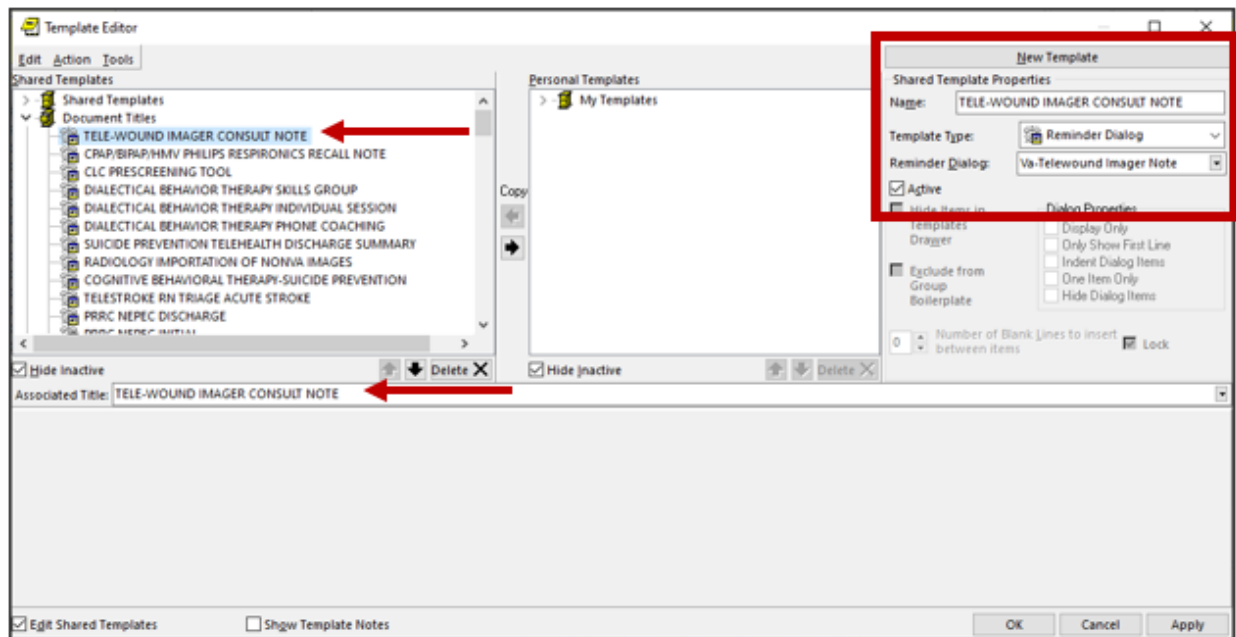


Type in the following information:

Name: **TELE-WOUND IMAGER CONSULT NOTE**

Template Type: **Reminder Dialog**

Reminder Dialog: **VA-TELEWOUND IMAGER NOTE**

Next, type in the Associated Title name – **TELE-WOUND IMAGER CONSULT NOTE** – then click Apply.

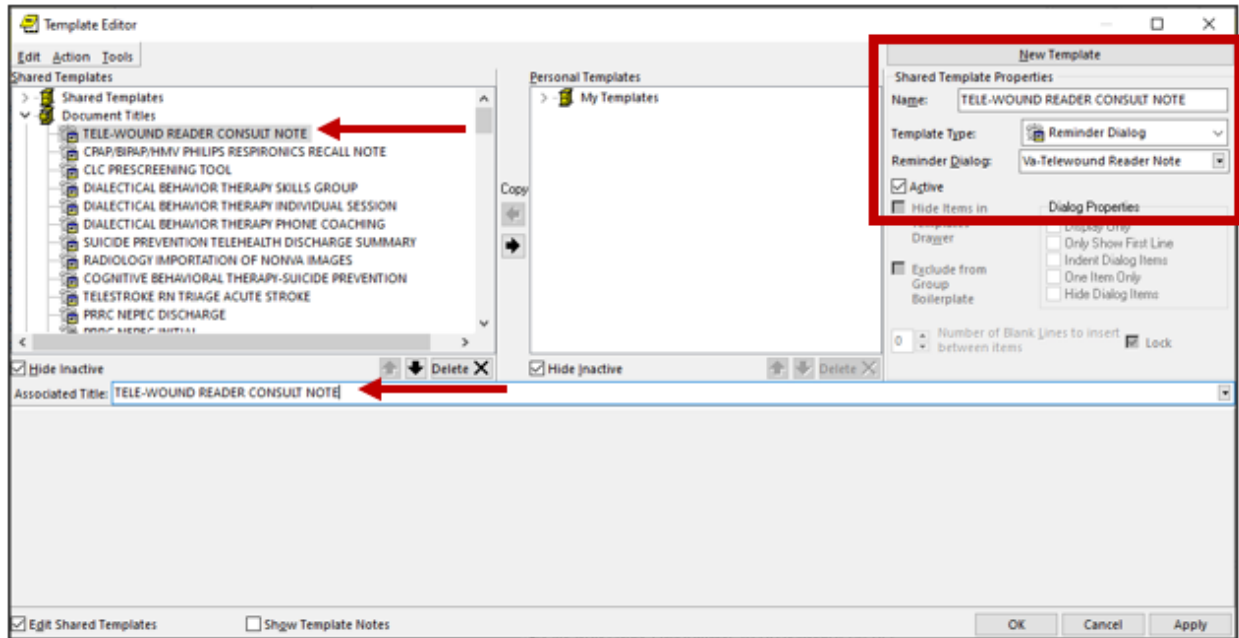

Type in the following information:

Name: **TELE-WOUND READER CONSULT NOTE**

Template Type: **Reminder Dialog**

Reminder Dialog: **VA-TELEWOUND READER NOTE**

Next, type in the Associated Title name – **TELE-WOUND READER CONSULT NOTE** – then click Apply.