



**Update \_2\_0\_316A**

**Clinical Reminders**

**VA-WH BREAST CANCER RISK ASSESSMENT  
Reminder Dialog Update**

**Install Guide**

**July 2023**

Product Development  
Office of Information Technology  
Department of Veterans Affairs



## **Table of Contents**

<b>INTRODUCTION .....</b>	<b>1</b>
<b>INSTALL DETAILS .....</b>	<b>2</b>
<b>INSTALL EXAMPLE .....</b>	<b>3</b>
<b>POST INSTALLATION .....</b>	<b>6</b>



# Introduction

---

## Description:

This update revises the VA-WH BREAST CANCER RISK ASSESSMENT reminder dialog released February 2023 in Update 316.

Installation is mandatory, but implementation and use are optional. Implementation should be coordinated with primary care and/or women's health staff.

Revisions: [316A Revisions](#)

Link to February 2023 Update 316 Education and Information: [Update 316](#)

UPDATE\_2\_0\_316A contains 1 Reminder Exchange entry:

## **UPDATE\_2\_0\_316A VA-BREAST CANCER RISK ASSESSMENT**

**The exchange file contains the following components:**

### **TIU TEMPLATE FIELD**

WORD PROCESSING (REQ)  
NCI BREAST CA RISK TOOL

### **HEALTH FACTORS**

VA-UPDATE\_2\_0\_316A  
VA-WH BREAST CA ASHKENAZI JEWISH ANCESTRY  
VA-WH BREAST CA FH RISK ASSMNT NEGATIVE  
VA-WH BREAST CA GENETIC COUNSELING REFER NO  
VA-WH BREAST CA GENETIC COUNSELING REFER YES  
VA-WH BREAST CA HX AGE<51 BIOLOGIC DAUGHTER  
VA-WH BREAST CA HX AGE<51 BIOLOGIC MATERNAL AUNT  
VA-WH BREAST CA HX AGE<51 BIOLOGIC MGM  
VA-WH BREAST CA HX AGE<51 BIOLOGIC MOTHER  
VA-WH BREAST CA HX AGE<51 BIOLOGIC PATERNAL AUNT  
VA-WH BREAST CA HX AGE<51 BIOLOGIC PGM  
VA-WH BREAST CA HX AGE<51 BIOLOGIC SISTER  
VA-WH BREAST CA HX AGE<51 PERSONAL  
VA-WH BREAST CA HX AGE>50 2 OR MORE CASES SAME SIDE FAM  
VA-WH BREAST CA HX BIOL MALE RELATIVE  
VA-WH BREAST CA PERSONAL/FH RISK ASSMNT TOTAL  
VA-WH BREAST CANCER RISK ASSESSMENT [C]  
VA-WH OVARIAN CA/FTC/PPC HX BIOLOGIC DAUGHTER  
VA-WH OVARIAN CA/FTC/PPC HX BIOLOGIC MATERNAL AUNT  
VA-WH OVARIAN CA/FTC/PPC HX BIOLOGIC MGM  
VA-WH OVARIAN CA/FTC/PPC HX BIOLOGIC MOTHER  
VA-WH OVARIAN CA/FTC/PPC HX BIOLOGIC PATERNAL AUNT

VA-WH OVARIAN CA/FTC/PPC HX BIOLOGIC PGM  
VA-WH OVARIAN CA/FTC/PPC HX BIOLOGIC SISTER  
VA-WH OVARIAN CA/FTC/PPC HX PERSONAL  
VA-WH PERSONAL HISTORY ADH  
VA-WH PERSONAL HX BREAST CANCER  
VA-WH PERSONAL HX CHEST WALL RADIATION 10-30YO  
VA-WH PERSONAL HX DCIS  
VA-WH PERSONAL HX LCIS  
VA-WH PERSONAL/FDR BREAST CA HX NEGATIVE  
VA-WH PERSONAL/FDR HX BRCA 1/BCRA 2  
VA-WH PERSONAL/FDR HX BRRS  
VA-WH PERSONAL/FDR HX COWDEN SYNDROME  
VA-WH PERSONAL/FDR HX LI-FRAUMENI SYNDROME  
VA-WH PERSONAL/FDR HX OTHER BC GENETIC CONDITIONS  
WH MAMMOGRAM [C]  
WH ORDER MAMMOGRAM SCREEN HF  
VA-WH BCRA SERVICE ACT AREA NO  
VA-WH BCRA TOXIC EXPOSURE POSITIVE  
VA-WH BCRA TOXIC EXPOSURE NEGATIVE  
VA-WH BCRA PERSONAL/FAMILY HISTORY POSITIVE  
VA-WH BCRA PERSONAL/FAMILY HISTORY NEGATIVE  
VA-WH PERSONAL/FDR BREAST CA HX NEGATIVE

### **REMINDER SPONSOR**

WOMEN VETERANS HEALTH PROGRAM

### **REMINDER TERM**

VA-REMINDER UPDATE\_2\_0\_316A  
VA-SEX MALE  
VA-WH BL AGE 30-34  
VA-AGE<30  
VA-WH BL AGE 35-39  
VA-WH BL TES REMINDER  
VA-WH BL AGE>34  
VA-WH BL TES REMINDER

### **REMINDER DIALOG**

VA-WH BREAST CANCER RISK ASSESSMENT

### **TIU OBJECTS**

VA-WH TES REMINDER RESULTS  
VA-WH BREAST CA RISK ASSMNT HX

## Pre-Installation

---

1. Identify local screening mammogram order or order menu. The order or menu will be mapped to the dialog during post-installation.
2. Identify local order for an interfacility consult to national TeleOncology (High Risk Breast Clinic), if available. The option in the dialog should be disabled if the interfacility consult is not available at your site. Instructions are in the post-installation steps.
3. Identify local order for consult to Women's Health. The order or menu will be mapped to the dialog during post-installation.

## Install Details

This update is being distributed as a web host file. The address for the host file is:

**https://REDACTED /UPDATE\_2\_0\_316A.PRD**

The file will be installed using Reminder Exchange, programmer access is not required.

Installation:

=====

This update can be loaded with users on the system. Installation will take less than **5** minutes.

## Install Example

---

To Load the Web Host File. Navigate to Reminder exchange in Vista

```
CFE Create Exchange File Entry      LHF Load Host File
CHF Create Host File                LMM Load MailMan Message
CMM Create MailMan Message          LR List Reminder Definitions
DFE Delete Exchange File Entry      LWH Load Web Host File
IFE Install Exchange File Entry      RI Reminder Definition Inquiry
IH Installation History              RP Repack
Select Action: Next Screen// LWH Load Web Host File
Input the URL for the .prd file: https://REDACTED /UPDATE_2_0_316A.PRD
```

At the **Select Action:** prompt, enter **LWH** for Load Web Host File

At the **Input the url for the .prd file:** prompt, type the following web address:

**https://REDACTED /UPDATE\_2\_0\_316A.PRD**

You should see a message at the top of your screen that the file successfully loaded.

```

497 UPDATE_2_0_316A VA-BREAST
    CANCER RISK ASSESSMENT
498 UPDATE_2_0_317 VA-CAREGIVER
    TEMPLATE UPDATES
499 UPDATE_2_0_318 VA-WHS
    WELL-BEING SIGNS UPDATE
500 UPDATE_2_0_319 VA-RPS VIRTUAL
    REALITY CREATIVE ARTS THERAPY
501 UPDATE_2_0_319 VA-RPS VIRTUAL
    REALITY CREATIVE ARTS THERAPY
+      + Next Screen  - Prev Screen  ?? More Actions  >>>
CFE Create Exchange File Entry      LHF Load Host File
CHF Create Host File                 LMM Load MailMan Message
CMM Create MailMan Message           LR  List Reminder Definitions
DFE Delete Exchange File Entry       LWH Load Web Host File
IFE Install Exchange File Entry      RI  Reminder Definition Inquiry
IH  Installation History              RP  Repack
Select Action: Next Screen// IFE    Install Exchange File Entry
Enter a list or range of numbers (1-798): 497

```

Search and locate an entry titled **UPDATE\_2\_0\_316A VA-BREAST CANCER RISK ASSESSMENT** in reminder exchange.

At the **Select Action** prompt, enter **IFE** for Install Exchange File Entry

Enter the number that corresponds with your entry UPDATE\_2\_0\_316A VA-BREAST CANCER RISK ASSESSMENT

(In this example it is entry 497. It will vary by site). The date of the exchange file should be 07/07/2023.

```

Description:
The following Clinical Reminder items were selected for packing:
REMINDER DIALOG
  VA-WH BREAST CANCER RISK ASSESSMENT

REMINDER TERM
  VA-REMINDER UPDATE_2_0_316A

Keywords:

Components:

+      Enter ?? for more actions  >>>
IA  Install all Components          IS  Install Selected Component
Select Action: Next Screen// ia    Install all Components  ....

```

At the **Select Action** prompt, type **IA** for Install all Components and hit enter.

Select Action: Next Screen// **IA Install all Components**

You will see several prompts, for all new entries, you will choose **I to Install**

**Overwrite** all items that already exist, but the packed component is different.



You will be prompted to install the reminder dialog component (VA-WH BREAST CANCER RISK ASSESSMENT). You will choose to install all even though some components may have an X to indicate that they exist.

```

Item  Seq.  Dialog Findings                                     Type  Exists
  1    VA-WH BREAST CANCER RISK ASSESSMENT                dialog  X
  2    5 VA-WH BCRA VERSION                               element
      Finding: *NONE*
  3    10 VA-WH BCRA SERVICE ACT AREA                    group
      Finding: *NONE*
  4    10.5 VA-WH BCRA SERVICE ACT AREA INFO              group
      Finding: VIEW PROGRESS NOTE TEXT (REMINDER GENERAL FINDING)  X
  5    PXRМ GF VIEW BUTTON                               prompt  X
  6    10.5.5 VA-WH BCRA SERVICE ACT AREA YES            group
      Finding: *NONE*
  7    10.5.5.2 VA-WH BCRA ENROLLED                      group
      Finding: *NONE*
  8    10.5.5.2.5 VA-WH BCRA ENROLLED YES                group
+      + Next Screen  - Prev Screen  ?? More Actions
DD  Dialog Details      DT  Dialog Text      IS  Install Selected
DF  Dialog Findings    DU  Dialog Usage    QU  Quit
DS  Dialog Summary     IA  Install All
Select Action: Next Screen// IA  Install All
Install reminder dialog and all components with no further changes: Y//

```

At the **Select Action** prompt, type **IA** to install the dialog: **VA-WH BREAST CANCER RISK ASSESSMENT**

Select Action: Next Screen// **IA Install All**

Install reminder dialog and all components with no further changes: Y// **Yes**

```

Packed reminder dialog: VA-WH BREAST CANCER RISK ASSESSMENT
VA-WH BREAST CANCER RISK ASSESSMENT (reminder dialog) installed from exchange fi
Item  Seq.  Dialog Findings                                     Type  Exists
  1    VA-WH BREAST CANCER RISK ASSESSMENT                dialog  X
  2    5 VA-WH BCRA VERSION                               element  X
      Finding: *NONE*
  3    10 VA-WH BCRA SERVICE ACT AREA                    group    X
      Finding: *NONE*
  4    10.5 VA-WH BCRA SERVICE ACT AREA INFO              group    X
      Finding: VIEW PROGRESS NOTE TEXT (REMINDER GENERAL FINDING)  X
  5    PXRМ GF VIEW BUTTON                               prompt   X
  6    10.5.5 VA-WH BCRA SERVICE ACT AREA YES            group    X
      Finding: *NONE*
  7    10.5.5.2 VA-WH BCRA ENROLLED                      group    X
      Finding: *NONE*
  8    10.5.5.2.5 VA-WH BCRA ENROLLED YES                group    X
+      + Next Screen  - Prev Screen  ?? More Actions
DF  Dialog Findings    DU  Dialog Usage    QU  Quit
DS  Dialog Summary     IA  Install All
Select Action: Next Screen// Q  Quit

```

After completing this dialog install, you will type **Q**.

```
Description:
The following Clinical Reminder items were selected for packing:
REMINDER DIALOG
  VA-WH BREAST CANCER RISK ASSESSMENT

REMINDER TERM
  VA-REMINDER UPDATE_2_0_316A

Keywords:

Components:

+      + Next Screen  - Prev Screen  ?? More Actions  >>>
IA  Install all Components          IS  Install Selected Component
Select Action: Next Screen// q  Q
```

You will then be returned to this screen. At the **Select Action** prompt, type **Q**.

Install complete.

# Post Installation

---

1. **Map the local order for an interfacility consult to national TeleOncology (High Risk Breast Clinic) to the following two elements, if available. Disable these elements if the interfacility consult is not available at your site.**
  - a. Navigate to the reminder dialog element menu in VistA and locate the elements.

NAME: VAL-WH BCRA PH RISK ASSMNT POSITIVE TELEONC  
Replace  
DISABLE: Enter "2" to disable this option if the consult to national TeleOncology (High Risk Breast Clinic) is not available.  
CLASS: LOCAL//  
SPONSOR:  
REVIEW DATE:  
RESOLUTION TYPE:  
ORDERABLE ITEM:  
FINDING ITEM: Enter "Q", then name of local quick order for the consult to national TeleOncology (High Risk Breast Clinic), if available.

NAME: VAL-WH BCRA FH RISK ASSMNT REFER TELEONC  
Replace  
DISABLE: Enter "2" to disable this option if the consult to national TeleOncology (High Risk Breast Clinic) is not available.  
CLASS: LOCAL//  
SPONSOR:  
REVIEW DATE:  
RESOLUTION TYPE:  
ORDERABLE ITEM:

FINDING ITEM: Enter "Q.", then name of local quick order for the consult to national TeleOncology (High Risk Breast Clinic), if available.

2. Map local screening mammogram order/menu to the VA-WH BCRA ORDER MAMMOGRAM element and VA-WH BCRA MAMMOGRAM ORDER DECISION element.

NAME: VA-WH BCRA ORDER MAMMOGRAM  
Finding item: HF WH ORDER MAMMOGRAM SCREEN HF <enter>  
FINDING ITEM: WH ORDER MAMMOGRAM SCREEN HF // <enter>

Additional findings: none  
Select ADDITIONAL FINDING: Enter "Q.", then name of local screening mammogram quick order or menu.

3. Map the local order for a consult to Women's Health for detailed risk assessment. If a consult doesn't exist, coordinate with Women's Health to change text in element to match local process.

- a. Navigate to the reminder dialog element menu in VistA and locate the elements.

NAME: VAL-WH BCRA PH RISK ASSMNT POSITIVE LOCAL REFER  
DISABLE:  
CLASS: LOCAL//  
SPONSOR:  
REVIEW DATE:  
RESOLUTION TYPE:  
ORDERABLE ITEM:  
FINDING ITEM: Enter "Q.", then name of local quick order for the consult to Women's Health.

Additional findings: none  
Select ADDITIONAL FINDING:  
DIALOG/PROGRESS NOTE TEXT:  
Positive - Referred for local detailed risk assessment

Edit? NO// YES Enter "Yes" if dialog text needs to be revised.

==[ WRAP ]==[INSERT ]===< DIALOG/PROGRESS NOTE TE[Press <PF1>H for help]====  
Positive - Referred for local detailed risk assessment

The text above can be revised, as needed.

<=====T=====T=====T=====T=====T=====T=====T=====T=====T=====T>=====

NAME: VAL-WH BCRA FH RISK ASSMNT REFER LOCAL  
DISABLE:  
CLASS: LOCAL//  
SPONSOR:

REVIEW DATE:  
RESOLUTION TYPE:  
ORDERABLE ITEM:  
FINDING ITEM: Enter "Q", then name of local quick order for the consult to Women's Health.

Additional findings: none  
Select ADDITIONAL FINDING:  
DIALOG/PROGRESS NOTE TEXT:  
Score of 2 or more- Referred for local/in person detailed risk assessment and genetic counseling

Edit? NO// YES Enter "Yes" if dialog text needs to be revised.

```
==[ WRAP ]==[INSERT ]===< DIALOG/PROGRESS NOTE TE[Press <PF1>H for help]====  
Score of 2 or more- Referred for local/in person detailed risk assessment  
and genetic counseling
```

The text above can be revised, as needed.

```
<=====T=====T=====T=====T=====T=====T=====T=====T=====T>=====
```

#### 4. Map the VA-TOXIC EXPOSURE reminder definition to the VA-WH BL TES REMINDER reminder term.

- b. Navigate to the Reminder Manager Menu in VistA. Follow the prompts below.

Select Reminder Dialog Management <TEST ACCOUNT> Option:

```
CF      Reminder Computed Finding Management ...  
RM      Reminder Definition Management ...  
SM      Reminder Sponsor Management ...  
TXM     Reminder Taxonomy Management  
TRM     Reminder Term Management ...  
LM      Reminder Location List Management ...  
RX      Reminder Exchange  
RT      Reminder Test  
OS      Other Supporting Menus ...  
INFO    Reminder Information Only Menu ...  
DM      Reminder Dialog Management ...  
CP      CPRS Reminder Configuration ...  
RP      Reminder Reports ...  
MST     Reminders MST Synchronization Management ...  
PL      Reminder Patient List Menu ...  
PAR     Reminder Parameters ...  
VS      NLM Value Set Menu  
ROC     Reminder Order Check Menu ...  
CQM     NLM Clinical Quality Measures Menu  
XM      Reminder Extract Menu ...
```

Select Reminder Managers Menu <TEST ACCOUNT> Option: TRM

```
TL      List Reminder Terms
```

TI Inquire about Reminder Term  
TE Add/Edit Reminder Term  
TC Copy Reminder Term  
TICS Integrity Check Selected  
TICA Integrity Check All  
TEST Term Test

Select Reminder Term Management <TEST ACCOUNT> Option: **TE**

Select Reminder **VA-WH BL TES REMINDER**  
...OK? Yes//**Yes**

NAME: VA-WH BL TES REMINDER **<enter>**  
CLASS: NATIONAL// **<enter>**  
SPONSOR: **<enter>**  
REVIEW DATE: **<enter>**  
DESCRIPTION: **<enter>**  
Edit? NO//

Choose from:

CF VA-REMINDER DEFINITION

Finding # 1

Select Finding: **\1**

**Below, Replace the existing reminder with the VA-TOXIC EXPOSURE reminder.**

Display help for CF.VA-REMINDER DEFINITION? **NO**

Editing Finding Number: 1 **<enter>**

FINDING ITEM: VA-REMINDER DEFINITION// **<enter>**

BEGINNING DATE/TIME: **<enter>**

ENDING DATE/TIME: **<enter>**

OCCURRENCE COUNT: **<enter>**

COMPUTED FINDING PARAMETER: VA-REPLACE DEFINITION

Replace ... With **VA-TOXIC EXPOSURE**

CONDITION: I V["DUE"**<enter>**

CONDITION CASE SENSITIVE: NO// **<enter>**

USE STATUS/COND IN SEARCH: **<enter>**

## 5. Embed the VA-TOXIC EXPOSURE dialog into the VAL-WH BCRA TOXIC EXPOSURE SCREEN group.

- a. Navigate to the reminder dialog group menu in VistA and locate the VAL-WH BCRA TOXIC EXPOSURE SCREEN group.

NAME: **VAL-WH BCRA TOXIC EXPOSURE SCREEN**  
DISABLE: **<enter>**  
CLASS: LOCAL/**<enter>**/  
SPONSOR: **<enter>**  
REVIEW DATE: **<enter>**  
RESOLUTION TYPE: **<enter>**  
ORDERABLE ITEM: **<enter>**  
FINDING ITEM: **<enter>**  
Additional findings: **<enter>**none  
Select ADDITIONAL FINDING: **<enter>**

GROUP CAPTION: <enter>  
 PUT A BOX AROUND THE GROUP: <enter>  
 SHARE COMMON PROMPTS: <enter>  
 MULTIPLE SELECTION: <enter>  
 HIDE/SHOW GROUP: HIDE// <enter>  
 GROUP HEADER DIALOG TEXT:  
 Environmental Military Exposure Screen\\  
 (use to document the Veteran's Military Exposure history)  
 Edit? NO// <enter>  
 GROUP HEADER ALTERNATE P/N TEXT:  
 Environmental Military Exposure Screen  
 Edit? NO// <enter>  
 EXCLUDE FROM PROGRESS NOTE: <enter>  
 SUPPRESS CHECKBOX: <enter>  
 NUMBER OF INDENTS: 2// <enter>  
 INDENT PROGRESS NOTE TEXT: INDENT//<enter>  
 Select SEQUENCE: 5  
 DIALOG ELEMENT: VA-TOXIC EXPOSURE SCREENING VERSION      dialog element  
     ...OK? Yes// Yes

EXCLUDE FROM PN TEXT: <enter>  
 Select SEQUENCE: 10  
 DIALOG ELEMENT: VAL-TOXIC EXPOSURE SCREEN NON CLINICAL      dialog element  
     ...OK? Yes// Yes

EXCLUDE FROM PN TEXT: <enter>  
 Select SEQUENCE: 15  
 DIALOG ELEMENT: VA-GP TOXIC EXPOSURE SCREEN  
     1 VA-GP TOXIC EXPOSURE SCREEN      dialog group  
     2 VA-GP TOXIC EXPOSURE SCREEN AHBP INFO      dialog element  
     3 VA-GP TOXIC EXPOSURE SCREEN ELIG QU CONTACTS      dialog group  
     4 VA-GP TOXIC EXPOSURE SCREEN ELIG QUESTIONS      dialog group  
     5 VA-GP TOXIC EXPOSURE SCREEN ENROLL QU CONTACTS      dialog group  
 CHOOSE 1-5: 1 VA-GP TOXIC EXPOSURE SCREEN      dialog group      NATIONAL  
 EXCLUDE FROM PN TEXT: <enter>

Select SEQUENCE: 15//<enter>  
 Select BRANCHING LOGIC SEQUENCE: 5// <enter>  
 BRANCHING LOGIC SEQUENCE: 5// <enter>  
 EVALUATION ITEM: VA-WH BL TES REMINDER// <enter>  
 EVALUATION STATUS: FALSE// <enter>  
 ACTION: REPLACE// <enter>  
 REPLACEMENT ITEM: VA-WH BCRA TES RESULTS//<enter>  
 Select BRANCHING LOGIC SEQUENCE: <enter>  
 Checking reminder dialog for errors..  
 NO ERRORS FOUND  
 Input your edit comments.  
 Edit? NO//<enter>

6. **If not completed after installing Update 316, make the dialog template available to be attached to a Personal or Shared template or to be able to attach the dialog to a progress note title.**

Select Reminder Managers Menu <TEST ACCOUNT> Option: **CP CPRS**

**Reminder Configuration**

- CA Add/Edit Reminder Categories
- CL CPRS Lookup Categories
- CS CPRS Cover Sheet Reminder List
- MH Mental Health Dialogs Active
- PN Progress Note Headers
- RA Reminder GUI Resolution Active
- TIU TIU Template Reminder Dialog Parameter**
- DEVL Evaluate Coversheet List on Dialog Finish
- DL Default Outside Location
- PT Position Reminder Text at Cursor
- LINK Link Reminder Dialog to Template
- TEST CPRS Coversheet Time Test
- NP New Reminder Parameters
- GEC GEC Status Check Active
- WH WH Print Now Active

Select CPRS Reminder Configuration <TEST ACCOUNT> Option: **TIU TIU**

**Template Reminder Dialog Parameter**

Reminder Dialogs allowed as Templates may be set for the following:

- |   |          |     |                               |
|---|----------|-----|-------------------------------|
| 1 | User     | USR | [choose from NEW PERSON]      |
| 3 | Service  | SRV | [choose from SERVICE/SECTION] |
| 4 | Division | DIV | [SALT LAKE CITY]              |
| 5 | System   | SYS | [NATREM.FO-SLC.MED.VA.GOV]    |

Enter selection: **5** *Choose the appropriate number for your site. Your site may do this by System or other levels. The example below uses SYSTEM level.*

Setting Reminder Dialogs allowed as Templates for System: NATREM.FO-SLC.MED.VA.GOV

Select Display Sequence: **?**

- 388 VA-COVID-19 TELEPHONE FOLLOW-UP
- 389 VA-SRA-CONSULT RESPONSE (D)

**When you type a question mark above, you will see the list of #'s (with dialogs) that are already taken. Choose a number **NOT** on this list. For this example, 390 is not present so I will use 390.**

Select Display Sequence: **390**

Are you adding 390 as a new Display Sequence? Yes// **YES**

Display Sequence: 390// <Enter> 390  
Clinical Reminder Dialog: **VA-WH BREAST CANCER RISK ASSESSMENT**  
reminder dialog NATIONAL  
...OK? Yes// <Enter> (Yes)

**7. If not completed after installing Update 316, set up note title.**

The note title to create to associate with this template is: **BREAST CANCER RISK ASSESSMENT**


VHA Enterprise Standard Title: **WOMENS HEALTH NOTE**

You can add this to the appropriate document class.

**Note:** Comprehensive information on Creating new Document Classes and Titles can be found beginning on page 54 of the TIU/ASU Implementation Guide at the following link: [http://www.va.gov/vdl/documents/Clinical/CPRS-Text\\_Integration\\_UTILITY\\_\(TIU\)/tiuim.pdf](http://www.va.gov/vdl/documents/Clinical/CPRS-Text_Integration_UTILITY_(TIU)/tiuim.pdf)

**8. If not completed after installing Update 316, associate the reminder dialogs with the note titles in CPRS**

In CPRS, open the Template Editor, expand Document Titles, select the appropriate folder (if applicable), then click New Template (top right-hand corner).

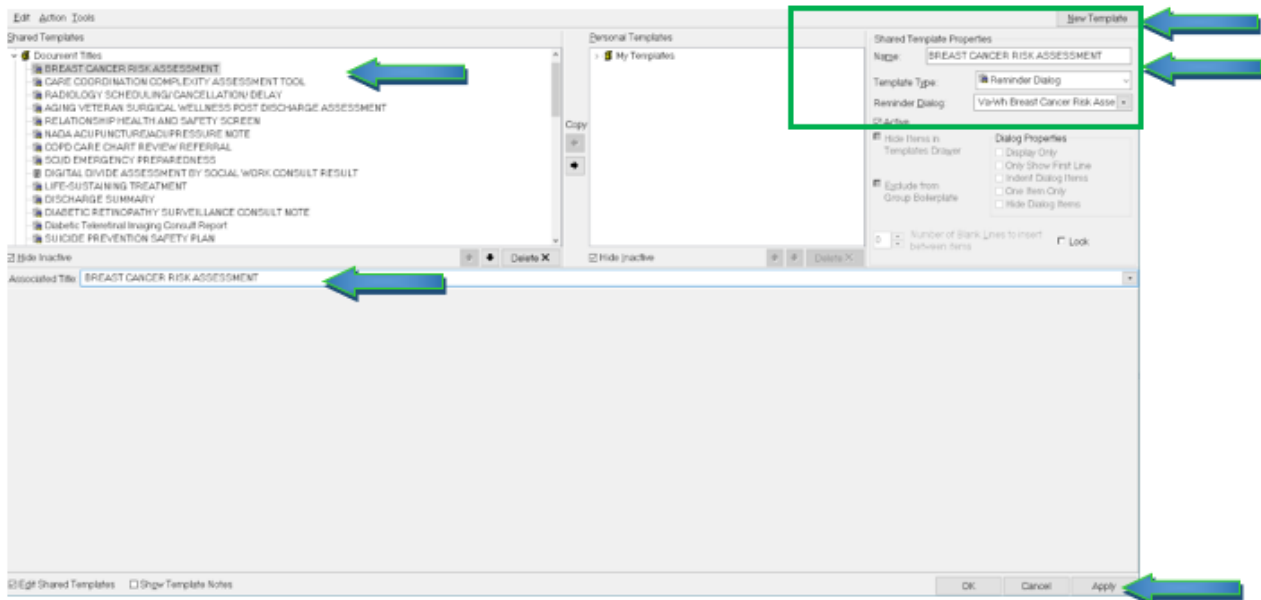

Type in the following information:



Name: **BREAST CANCER RISK ASSESSMENT**

Template Type: **Reminder Dialog**

Reminder Dialog: **VA-WH BREAST CANCER RISK ASSESSMENT**

Next, type in the Associated Title name — then click Apply.

9. **Associate the reminder dialog with a shared template in CPRS, if needed.**