



Update _2_0_370

Clinical Reminders

VA-TELERETINOPATHY Reminder and Dialog Updates

Install Guide

October 2023

Product Development
Office of Information Technology
Department of Veterans Affairs

Table of Contents

INTRODUCTION 1
INSTALL DETAILS 8
PRE INSTALLATION 9
INSTALL EXAMPLE 9
POST INSTALLATION 17

Introduction

DESCRIPTION:

This update revises several items in the VA-TELERETINOPATHY READER NOTE, VA-TELERETINOPATHY IMAGER NOTE, and VA-TELE-EYE RESULTS LETTER dialogs, the VA-EYE CARE AT-RISK SCREEN (NURSE), the VA-GP TECS TECH-DIABETES dialog group, and VA-EYE CARE AT-RISK SCREEN (PROV) reminder definitions. A detailed description of the revisions is listed below.

For a description of the reminder dialogs and definitions, see the original installation guides on the VA Software Documentation Library (VDL)

- UPDATE_2_0_167 Tele-Eye and Tele-Derm Updates
- UPDATE_2_0_146 VA-EYE CARE AT-RISK SCREEN
- UPDATE_2_0_42 Teleretinal Imaging Dialogs
- UPDATE_2_0_217 VA-Teleretinopathy Reminder Updates

*Note: The VA-Eye Care At-Risk Screen Reminder Definition may only be activated at sites prepared for expansion. Coordinate implementation on CPRS cover sheet with the Facility Telehealth Coordinator.

*Note: Use of the VA-EYE CARE AT-RISK SCREEN 'NURSE' reminder is by Primary Care Nursing and/or Telehealth Imager/staff.

*Note: Reminders may be installed and activated on a community-based outpatient clinic (CBOC) by CBOC basis.

Updates to the two reminders and five dialogs are listed here:

UPDATE_2_0_370 contains 1 Reminder Exchange entry:

UPDATE_2_0_370 VA-TELERETINOPATHY UPDATE

The exchange file contains the following components:

TEMPLATE FIELDS

VA-TECS DURATION OF DIABETES

WP FREE TEXT FL55 I3

COMMENT OPTIONAL

VA-TECS OD/OS/OU

BLANK SPACE1

VA-YES NO REQUIRED

DATE

VA-TRR DATE 9M

VA-TRR DATE 2Y

VA-TRR DATE 1Y

VA-TRR DATE 6M

VA-TRR DATE 4M

VA-TRR DATE 3M

VA-TRR DATE 2M
VA-TRR DATE 1M
TEXT (1-60 CHAR)

EXAM

TRR EXAM OTH ASSESS-U/A
TRR EXAM OTH ASSESS-ABN
TRR EXAM OTH ASSESS-NORM
TRR EXAM OPTC NRV-NORMAL
TRR EXAM OPTC NRV-ABN
TRR EXAM OPTC NRV-U/A
TRR EXAM MACULA-U/A
TRR EXAM MACULA-ABN
TRR EXAM MACULA-NORMAL
VA-TRR DM EXAM RTNOPTY-U/A
VA-TRR DM EXAM RTNOPTY-ABN
VA-TRR DM EXAM RTNOPTY-NORMAL

HEALTH FACTORS

INCORRECT DIAGNOSES [C]
LIFE EXPECTANCY < 6 MONTHS
REMINDER FACTORS [C]
VA-UPDATE_2_0_370
VA-COMP EYE EXAM SCHEDULED
VA-EYE AT RISK SCREEN HX LASER TX/INJ LAST 2YRS
VA-EYE CARE AT-RISK SCREEN [C]
VA-INCORRECT DIABETES DIAGNOSIS
VA-INCORRECT DIABETES DIAGNOSIS-PROVIDER
VA-OUTSIDE DIAB EYE ABN AT-RISK SCRIN
VA-OUTSIDE DIAB EYE NORMAL AT-RISK SCRIN
VA-OUTSIDE DIAB EYE UNK AT-RISK SCRIN
VA-OUTSIDE EYE AT-RISK SCREEN
VA-REFUSED NURS AT-RISK SCREEN
VA-REFUSED PROV AT-RISK SCREEN
VA-TECS EYE DX HX DM-INSULIN
VA-TECS EYE DX HX-DIABETES
VA-TECS TECHNICIAN [C]
VA-TELE EYE GLAUC EXAM UNABLE TO ASSESS
VA-TELE EYE GLAUCOMA EXAM ABNORMAL
VA-TELE EYE GLAUCOMA EXAM NORMAL
VA-TRI CATEGORY-DM EYE EXAM [C]
VA-TRI DM RETINOP EXAM ABNORMAL
VA-TRI DM RETINOP EXAM NORMAL
VA-TRI DM RETINOP EXAM UNABLE TO ASSESS
VA-TRI EYE EXAM FUTURE
VA-TRI EYE EXAM SCHEDULED
VA-TRI MACULA EXAM ABNORMAL
VA-TRI MACULA EXAM NORMAL
VA-TRI MACULA EXAM UNABLE TO ASSESS
VA-TRR AT RISK FOR GLAUCOMA
VA-TRR AT RISK FOR MAC DEGEN & GLAUCOMA
VA-TRR AT RISK FOR MAC DEGENERATION
VA-TRR AT RISK OTHER EYE CONDITION
VA-TRR CATEGORY-DM EXAM FINDINGS LT EYE [C]
VA-TRR CATEGORY-DM EXAM FINDINGS RT EYE [C]
VA-TRR CATEGORY-DM EYE EXAM [C]
VA-TRR CATEGORY-EYE SCREENING EXAM FINDINGS [C]
VA-TRR CATEGORY-IMAGE QUALITY [C]

VA-TRR CATEGORY-NON DM EYE EXAM [C]
VA-TRR CATEGORY-NON DM FINDINGS LT EYE [C]
VA-TRR CATEGORY-NON DM FINDINGS RT EYE [C]
VA-TRR CATEGORY-RECOMMENDATIONS [C]
VA-TRR CATEGORY-RESULTS NOTIFICATION [C]
VA-TRR CATEGORY-SCREEN [C]
VA-TRR DIABETES TYPE 1
VA-TRR DIABETES TYPE 2
VA-TRR DM MACULA EXAM UNABLE TO ASSESS
VA-TRR DM OPTIC NRV EXAM NORMAL
VA-TRR DM RETINOP EXAM ABNORMAL
VA-TRR DM RETINOP EXAM NORMAL
VA-TRR DM RETINOP EXAM UNABLE TO ASSESS
VA-TRR DM1 FIND DR MLD MAC NO EDEMA LT
VA-TRR DM1 FIND DR MLD MAC NO EDEMA RT
VA-TRR DM1 FIND DR MOD MAC NO EDEMA LT
VA-TRR DM1 FIND DR MOD MAC NO EDEMA RT
VA-TRR DM1 FIND DR PDR MAC NO EDEMA LT
VA-TRR DM1 FIND DR PDR MAC NO EDEMA RT
VA-TRR DM1 FIND DR PREV RET LASER TX LT
VA-TRR DM1 FIND DR PREV RET LASER TX RT
VA-TRR DM1 FIND DR SEV MAC NO EDEMA LT
VA-TRR DM1 FIND DR SEV MAC NO EDEMA RT
VA-TRR DM1 FINDING DR INADQ ASSESS LT
VA-TRR DM1 FINDING DR INADQ ASSESS RT
VA-TRR DM1 FINDING DR MLD MAC EDEMA LT
VA-TRR DM1 FINDING DR MLD MAC EDEMA RT
VA-TRR DM1 FINDING DR MOD MAC EDEMA LT
VA-TRR DM1 FINDING DR MOD MAC EDEMA RT
VA-TRR DM1 FINDING DR NO RETINOPATHY LT
VA-TRR DM1 FINDING DR NO RETINOPATHY RT
VA-TRR DM1 FINDING DR OTHER RT
VA-TRR DM1 FINDING DR PDR MAC EDEMA LT
VA-TRR DM1 FINDING DR PDR MAC EDEMA RT
VA-TRR DM1 FINDING DR SEV MAC EDEMA LT
VA-TRR DM1 FINDING DR SEV MAC EDEMA RT
VA-TRR DM2 FIND DR MLD MAC NO EDEMA LT
VA-TRR DM2 FIND DR MLD MAC NO EDEMA RT
VA-TRR DM2 FIND DR MOD MAC NO EDEMA LT
VA-TRR DM2 FIND DR MOD MAC NO EDEMA RT
VA-TRR DM2 FIND DR PDR MAC NO EDEMA LT
VA-TRR DM2 FIND DR PDR MAC NO EDEMA RT
VA-TRR DM2 FIND DR PREV RET LASER TX LT
VA-TRR DM2 FIND DR PREV RET LASER TX RT
VA-TRR DM2 FIND DR SEV MAC NO EDEMA LT
VA-TRR DM2 FIND DR SEV MAC NO EDEMA RT
VA-TRR DM2 FINDING DR INADQ ASSESS LT
VA-TRR DM2 FINDING DR INADQ ASSESS RT
VA-TRR DM2 FINDING DR MLD MAC EDEMA LT
VA-TRR DM2 FINDING DR MLD MAC EDEMA RT
VA-TRR DM2 FINDING DR MOD MAC EDEMA LT
VA-TRR DM2 FINDING DR MOD MAC EDEMA RT
VA-TRR DM2 FINDING DR NO RETINOPATHY LT
VA-TRR DM2 FINDING DR NO RETINOPATHY RT
VA-TRR DM2 FINDING DR OTHER LT
VA-TRR DM2 FINDING DR OTHER RT
VA-TRR DM2 FINDING DR PDR MAC EDEMA LT
VA-TRR DM2 FINDING DR PDR MAC EDEMA RT
VA-TRR DM2 FINDING DR SEV MAC EDEMA LT

VA-TRR DM2 FINDING DR SEV MAC EDEMA RT
VA-TRR IMAGE QUAL ADEQUATE
VA-TRR IMAGE QUAL INADEQUATE
VA-TRR IMAGE QUAL INADQ-BLURRED IMAGE
VA-TRR IMAGE QUAL INADQ-FIELD MISSING
VA-TRR IMAGE QUAL INADQ-MEDIA
VA-TRR IMAGE QUAL INADQ-OTHER REASON
VA-TRR IMAGE QUAL INADQ-PT POSITION
VA-TRR IMAGE QUAL INADQ-PUPIL SIZE
VA-TRR MAC FINDING ATROPHY LT
VA-TRR MAC FINDING ATROPHY RT
VA-TRR MAC FINDING DRUSEN LT
VA-TRR MAC FINDING DRUSEN RT
VA-TRR MAC FINDING DRUSEN-LG LT
VA-TRR MAC FINDING DRUSEN-LG RT
VA-TRR MAC FINDING DRUSEN-SMALL LT
VA-TRR MAC FINDING DRUSEN-SMALL RT
VA-TRR MAC FINDING EXUDATES LT
VA-TRR MAC FINDING EXUDATES RT
VA-TRR MAC FINDING GEOG ATROPHY LT
VA-TRR MAC FINDING GEOG ATROPHY RT
VA-TRR MAC FINDING INADEQ ASSESS LT
VA-TRR MAC FINDING INADEQ ASSESS RT
VA-TRR MAC FINDING INTRARETINAL HEMORRHAGE LT
VA-TRR MAC FINDING INTRARETINAL HEMORRHAGE RT
VA-TRR MAC FINDING MAC HOLE-LT
VA-TRR MAC FINDING MAC HOLE-RT
VA-TRR MAC FINDING NEOVASC MACULOPHTHY LT
VA-TRR MAC FINDING NEOVASC MACULOPHTHY RT
VA-TRR MAC FINDING NO APPARENT ABN LT
VA-TRR MAC FINDING NO APPARENT ABN RT
VA-TRR MAC FINDING OTHER LT
VA-TRR MAC FINDING OTHER RT
VA-TRR MAC FINDING PIGMENTARY CHANGE LT
VA-TRR MAC FINDING PIGMENTARY CHANGE RT
VA-TRR MAC FINDING PUCKERING LT
VA-TRR MAC FINDING PUCKERING RT
VA-TRR MACULA EXAM ABNORMAL
VA-TRR MACULA EXAM NORMAL
VA-TRR MACULA EXAM UNABLE TO ASSESS
VA-TRR NOT CANDIDATE FOR SCRIN
VA-TRR ON FINDING DISC DRUSEN LT
VA-TRR ON FINDING DISC DRUSEN RT
VA-TRR ON FINDING DISC EDEMA SUSP LT
VA-TRR ON FINDING DISC EDEMA SUSP RT
VA-TRR ON FINDING DISC NEOVASCULARIZATION LT
VA-TRR ON FINDING DISC NEOVASCULARIZATION RT
VA-TRR ON FINDING DISC W/ SHARP MARGINS LT
VA-TRR ON FINDING DISC W/ SHARP MARGINS RT
VA-TRR ON FINDING ELEVATED IOP LT
VA-TRR ON FINDING ELEVATED IOP RT
VA-TRR ON FINDING INADQ ASSESS LT
VA-TRR ON FINDING INADQ ASSESS RT
VA-TRR ON FINDING KNOWN GLAUCOMA LT
VA-TRR ON FINDING KNOWN GLAUCOMA RT
VA-TRR ON FINDING LG CUP/DISC LT
VA-TRR ON FINDING LG CUP/DISC RT
VA-TRR ON FINDING NO APP ABN LT
VA-TRR ON FINDING NO APP ABN RT
VA-TRR ON FINDING OTHER LT
VA-TRR ON FINDING OTHER RT

VA-TRR ON FINDING PALLOR LT
VA-TRR ON FINDING PALLOR RT
VA-TRR ON FINDING PAPILLEDEMA SUSPECTED LT
VA-TRR ON FINDING PAPILLEDEMA SUSPECTED RT
VA-TRR ON FINDING PERIPAPILLARY ATROPHY LT
VA-TRR ON FINDING PERIPAPILLARY ATROPHY RT
VA-TRR ON FINDING PHYSIOLOGIC CUPPING SUSPECTED LT
VA-TRR ON FINDING PHYSIOLOGIC CUPPING SUSPECTED RT
VA-TRR OPTIC NRV EXAM ABNORMAL
VA-TRR OPTIC NRV EXAM NORMAL
VA-TRR OPTIC NRV EXAM UNABLE TO ASSES
VA-TRR OTH ASSESS EXAM ABNORMAL
VA-TRR OTH ASSESS EXAM NORMAL
VA-TRR OTH ASSESS EXAM UNABLE TO ASSESS
VA-TRR OTH ASSESS FIND INADQ ASSESS LT
VA-TRR OTH ASSESS FIND INADQ ASSESS RT
VA-TRR OTH ASSESS FINDING CATARACT LT
VA-TRR OTH ASSESS FINDING CATARACT RT
VA-TRR OTH ASSESS FINDING CHOR NEVUS LT
VA-TRR OTH ASSESS FINDING CHOR NEVUS RT
VA-TRR OTH ASSESS FINDING CHORIORETINAL SCAR LT
VA-TRR OTH ASSESS FINDING CHORIORETINAL SCAR RT
VA-TRR OTH ASSESS FINDING COTTON WOOL SPOT LT
VA-TRR OTH ASSESS FINDING COTTON WOOL SPOT RT
VA-TRR OTH ASSESS FINDING FAMILY HX GLAUCOMA LT
VA-TRR OTH ASSESS FINDING FAMILY HX GLAUCOMA RT
VA-TRR OTH ASSESS FINDING HTN RTNPHTY LT
VA-TRR OTH ASSESS FINDING HTN RTNPHTY RT
VA-TRR OTH ASSESS FINDING LENS IMPLANT LT
VA-TRR OTH ASSESS FINDING LENS IMPLANT RT
VA-TRR OTH ASSESS FINDING LENTICULAR OPACIFICATION LT
VA-TRR OTH ASSESS FINDING LENTICULAR OPACIFICATION RT
VA-TRR OTH ASSESS FINDING NO ABN FIND LT
VA-TRR OTH ASSESS FINDING NO ABN FIND RT
VA-TRR OTH ASSESS FINDING OCCLUSION LT
VA-TRR OTH ASSESS FINDING OCCLUSION RT
VA-TRR OTH ASSESS FINDING OTH LT
VA-TRR OTH ASSESS FINDING OTH RT
VA-TRR OTH ASSESS FINDING PERIPH RET DEGEN LT
VA-TRR OTH ASSESS FINDING PERIPH RET DEGEN RT
VA-TRR OTH ASSESS FINDING RET EMBOLUS LT
VA-TRR OTH ASSESS FINDING RET EMBOLUS RT
VA-TRR OTH ASSESS FINDING RETINAL VASC CHANGE LT
VA-TRR OTH ASSESS FINDING RETINAL VASC CHANGE RT
VA-TRR RECOMMEND VA EYE CLINIC APPT
VA-TRR REFER COMP EYE CARE
VA-TRR REFER MACULA DISEASE
VA-TRR REFER OPHTH/OPTOM 1 MONTH
VA-TRR REFER OPHTH/OPTOM 1 YR
VA-TRR REFER OPHTH/OPTOM 2 MONTHS
VA-TRR REFER OPHTH/OPTOM 2 YRS
VA-TRR REFER OPHTH/OPTOM 3 MONTHS
VA-TRR REFER OPHTH/OPTOM 4 MONTHS
VA-TRR REFER OPHTH/OPTOM 6 MONTHS
VA-TRR REFER OPHTH/OPTOM 9 MONTHS
VA-TRR REFER OPHTH/OPTOM ASAP
VA-TRR REFER OPHTH/OPTOM-OTHER FREQ
VA-TRR REFER OPHTHALMOLOGY 1 MONTH
VA-TRR REFER OPHTHALMOLOGY 1 YR

VA-TRR REFER OPHTHALMOLOGY 2 MONTHS
VA-TRR REFER OPHTHALMOLOGY 2 YRS
VA-TRR REFER OPHTHALMOLOGY 3 MONTHS
VA-TRR REFER OPHTHALMOLOGY 4 MONTHS
VA-TRR REFER OPHTHALMOLOGY 6 MONTHS
VA-TRR REFER OPHTHALMOLOGY 9 MONTHS
VA-TRR REFER OPHTHALMOLOGY ASAP
VA-TRR REFER OPHTHALMOLOGY-OTHER FREQ
VA-TRR REFER OPTIC NERVE DISEASE
VA-TRR REFER OPTOMETRY 1 MONTH
VA-TRR REFER OPTOMETRY 1 YR
VA-TRR REFER OPTOMETRY 2 MONTHS
VA-TRR REFER OPTOMETRY 2 YRS
VA-TRR REFER OPTOMETRY 3 MONTHS
VA-TRR REFER OPTOMETRY 4 MONTHS
VA-TRR REFER OPTOMETRY 6 MONTHS
VA-TRR REFER OPTOMETRY 9 MONTHS
VA-TRR REFER OPTOMETRY ASAP
VA-TRR REFER OPTOMETRY-OTHER FREQ
VA-TRR REFER OTHER EYE DISEASE
VA-TRR REFER OUTSIDE PROV 1 MONTH
VA-TRR REFER OUTSIDE PROV 1 YR
VA-TRR REFER OUTSIDE PROV 2 MONTHS
VA-TRR REFER OUTSIDE PROV 2 YRS
VA-TRR REFER OUTSIDE PROV 3 MONTHS
VA-TRR REFER OUTSIDE PROV 4 MONTHS
VA-TRR REFER OUTSIDE PROV 6 MONTHS
VA-TRR REFER OUTSIDE PROV 9 MONTHS
VA-TRR REFER OUTSIDE PROV ASAP
VA-TRR REFER RETINOPATHY
VA-TRR REFER UNABLE TO ASSESS
VA-TRR REFER VISUALLY SIGNF CATARACT
VA-TRR REPEAT IMAGING-RET NEXT AVAIL
VA-TRR REPEAT IMAGING-RETURN 1 YR
VA-TRR REPEAT IMAGING-RETURN 2 YRS
VA-TRR REPEAT IMAGING-RETURN OTH
VA-TRR RESULTS NOTIFICATION COMPLETED

REMINDER SPONSOR

OFFICE OF TELEHEALTH SERVICES STORE AND FORWARD TELEHEALTH
OFFICE OF PATIENT CARE SERVICES
NATIONAL CENTER FOR HEALTH PROMOTION AND DISEASE PREVENTION (NCP)

REMINDER TERM

TELE EYE GLAUCOMA EXAM-ABNORMAL
TELE EYE GLAUCOMA EXAM-NORMAL
TELE EYE GLAUCOMA EXAM-UNABLE TO ASSESS
TELE EYE MACULA EXAM-ABNORMAL
TELE EYE MACULA EXAM-NORMAL
TELE EYE MACULA EXAM-UNABLE TO ASSESS
TELE EYE RETINOPATHY EXAM-ABNORMAL
TELE EYE RETINOPATHY EXAM-NORMAL
TELE EYE RETINOPATHY EXAM-UNABLE TO ASSESS
TELE EYE SCREEN ORDERS
VA-AMD RISK FACTORS
VA-BL TECS DIABETES HX
VA-BL TELE EYE AMD PT
VA-BL TELE EYE DIABETIC PT
VA-BL TELE EYE GLAUCOMA PT
VA-BL TELE EYE GLAUCOMA/AMD PT

VA-BL TELE EYE INCORRECT DIABETIC DX
VA-EYE AT RISK NSG SCREEN INCORRECT DM DIAGNOSIS
VA-EYE AT RISK NSG SCREEN REFUSAL
VA-GLAUCOMA RISK FACTORS
VA-LIFE EXPECTANCY < 6 MONTHS
VA-OUTSIDE EYE AT-RISK SCREEN
VA-TRI SCREEN EXCLUSIONS
VA-REMINDER UPDATE_2_0_370

REMINDER TAXONOMY

LONG TERM STEROID USE
TOBACCO USE CODES
TRI EX DM RETINOPATHY U/A
VA-ABNORMAL BMI
VA-AMD RISK FACTORS
VA-CAD
VA-DIABETES
VA-DIABETES HEDIS
VA-DIABETIC W/ EYE CONDITIONS
VA-DIABETIC W/O EYE CONDITIONS
VA-EYE AT-RISK SCR N CPT EXCLUSIONS
VA-EYE AT-RISK SCR N ICD EXCLUSIONS
VA-GLAUCOMA RISK FACTORS
VA-GLAUCOMA/AMD W EYE CONDITION
VA-GLAUCOMA/AMD W/O EYE CONDITION
VA-TERMINAL CANCER PATIENTS
VA-TRI BRANCH RETINAL VEIN OCCLUSION
VA-TRI CHOROIDAL NEVUS L
VA-TRI CHOROIDAL NEVUS R
VA-TRI DM1 MILD NONPROLIF RETINOPATHY NO MACULAR EDEMA
VA-TRI DM1 MILD NONPROLIF RETINOPATHY NO MACULAR EDEMA R
VA-TRI DM1 MILD NONPROLIF RETINOPATHY PROB MACULAR EDEMA L
VA-TRI DM1 MILD NONPROLIF RETINOPATHY PROB MACULAR EDEMA R
VA-TRI DM1 MOD NONPROLIF RETINOPATHY NO MACULAR EDEMA L
VA-TRI DM1 MOD NONPROLIF RETINOPATHY NO MACULAR EDEMA R
VA-TRI DM1 MOD NONPROLIF RETINOPATHY PROB MACULAR EDEMA L
VA-TRI DM1 MOD NONPROLIF RETINOPATHY PROB MACULAR EDEMA R
VA-TRI DM1 PROLIF RETINOPATHY NO MACULAR EDEMA L
VA-TRI DM1 PROLIF RETINOPATHY NO MACULAR EDEMA R
VA-TRI DM1 PROLIF RETINOPATHY PROB MACULAR EDEMA L
VA-TRI DM1 PROLIF RETINOPATHY PROB MACULAR EDEMA R
VA-TRI DM1 SEVERE NONPROLIF RETINOPATHY NO MACULAR EDEMA L
VA-TRI DM1 SEVERE NONPROLIF RETINOPATHY NO MACULAR EDEMA R
VA-TRI DM1 SEVERE NONPROLIF RETINOPATHY PROB MACULAR EDEMA L
VA-TRI DM1 SEVERE NONPROLIF RETINOPATHY PROB MACULAR EDEMA R
VA-TRI DM1/DM2 NO RETINOPATHY NO MACULAR EDEMA R/L
VA-TRI DM2 MILD NONPROLIF RETINOPATHY NO MACULAR EDEMA L
VA-TRI DM2 MILD NONPROLIF RETINOPATHY NO MACULAR EDEMA R
VA-TRI DM2 MILD NONPROLIF RETINOPATHY PROB MACULAR EDEMA L
VA-TRI DM2 MILD NONPROLIF RETINOPATHY PROB MACULAR EDEMA R
VA-TRI DM2 MOD NONPROLIF RETINOPATHY NO MACULAR EDEMA L
VA-TRI DM2 MOD NONPROLIF RETINOPATHY NO MACULAR EDEMA R
VA-TRI DM2 MOD NONPROLIF RETINOPATHY PROB MACULAR EDEMA L
VA-TRI DM2 MOD NONPROLIF RETINOPATHY PROB MACULAR EDEMA R
VA-TRI DM2 PROLIF RETINOPATHY NO MACULAR EDEMA L
VA-TRI DM2 PROLIF RETINOPATHY NO MACULAR EDEMA R
VA-TRI DM2 PROLIF RETINOPATHY PROB MACULAR EDEMA L
VA-TRI DM2 PROLIF RETINOPATHY PROB MACULAR EDEMA R

VA-TRI DM2 SEVERE NONPROLIF RETINOPATHY NO MACULAR EDEMA L
VA-TRI DM2 SEVERE NONPROLIF RETINOPATHY NO MACULAR EDEMA R
VA-TRI DM2 SEVERE NONPROLIF RETINOPATHY PROB MACULAR EDEMA L
VA-TRI DM2 SEVERE NONPROLIF RETINOPATHY PROB MACULAR EDEMA R
VA-TRI ELEVATED IOP L

REMINDER DEFINITION

VA-OBJ AMD RISK FACTORS
VA-OBJ GLAUCOMA RISK FACTORS
VA-DIABETES OBJECT
VA-TRI SCREENING EXCLUSIONS
VA-BL TELE-EYE SCREEN-AMD
VA-BL TELE-EYE SCREEN-GLAUCOMA
VA-BL TELE-EYE SCREEN-GLAUCOMA/AMD RISKS
VA-BL TELE-EYE SCREEN-INCORRECT DM DX
VA-EYE CARE AT-RISK SCREEN (PROV)
VA-EYE CARE AT-RISK SCREEN (NURSE)

REMINDER DIALOG

VA-EYE CARE AT-RISK SCREEN (NURSING)
VA-EYE CARE AT-RISK SCREEN (PROV)
VA-TELERETINOPATHY IMAGER NOTE
VA-TELERETINOPATHY READER NOTE
VA-TELE-EYE RESULTS LETTER

HEALTH SUMMARY TYPE

VA-TECS DIABETES HX
VA-TELE EYE RECOMMENDATIONS
VA-TELE EYE SCREENING RESULTS
VA-AMD RISK FACTORS
VA-GLAUCOMA RISK FACTORS
VA-DIABETES OBJECT
VA-TRI SCREEN EXCLUSIONS
VA-TRI HGBA1C

TIU DOCUMENT DEFINITION

VA-TECS DIABETES HX
VA-TELE EYE RECOMMENDATIONS
VA-TELE EYE SCREENING RESULTS
VA-AMD RISK FACTORS
VA-GLAUCOMA RISK FACTORS
VA-DIABETES OBJECT
VA-TRI SCRIN EXCLUSIONS

Install Details

This update is being distributed as a web host file. The address for the host file is:

https://REDACTED/UPDATE_2_0_370.PRD

The file will be installed using Reminder Exchange, programmer access is not required.

Installation:

=====

This update can be loaded with users on the system. Installation will take less than 30 minutes.

Pre Installation

1. **Print/Screen capture the following terms to capture any mapped local findings. Terms should be merged during installation. The copy will be helpful to confirm the findings were maintained post installation.**
 - a. TELE EYE GLAUCOMA EXAM-ABNORMAL
 - b. TELE EYE GLAUCOMA EXAM-NORMAL
 - c. TELE EYE GLAUCOMA EXAM-UNABLE TO ASSESS
 - d. TELE EYE MACULA EXAM-ABNORMAL
 - e. TELE EYE MACULA EXAM-NORMAL
 - f. TELE EYE MACULA EXAM-UNABLE TO ASSESS
 - g. TELE EYE RETINOPATHY EXAM-ABNORMAL
 - h. TELE EYE RETINOPATHY EXAM-NORMAL
 - i. TELE EYE RETINOPATHY EXAM-UNABLE TO ASSESS
 - j. TELE EYE SCREEN ORDERS
 - k. VA-OUTSIDE EYE AT-RISK SCREEN
 - l. VA-EYE AT RISK NSG SCREEN REFUSAL
 - m. VA-EYE AT RISK NSG SCREEN INCORRECT DM DIAGNOSIS
 - n. VA-BL TELE EYE INCORRECT DIABETIC DX
2. **Identify orders mapped to the following elements as additional findings. The orders will need to be re-mapped after installation.**
 - a. Element VAL-OI EYE CARE AT-RISK RETINAL IMAGE CONSULT
 - b. Element VAL-OI EYE CARE AT-RISK EYE GLASS CONSULT
 - c. Element VAL-OI EYE CARE AT-RISK COMP EYE EXAM CONSULT
3. **Identify local consult order for nutrition services. It will be mapped to an element in post-installation.**
4. **Identify local TIU object that displays the most recent eye exam dates, if available.**
5. **Print/Screen capture the Health Summary Type TRI HGBA1C to capture mapped local hemoglobin A1C lab tests. The copy will be helpful to confirm the lab tests were maintained post installation.**
6. **Note: The VA-Eye Care At-Risk Screen Reminder Definitions were only activated at expansion locations after the Facility Telehealth Coordinator (FTC) notified the CAC when to activate. Coordinate activation of these cover sheet reminders with the local FTC, if not already active.**

Install Example

To Load the Web Host File. Navigate to Reminder exchange in Vista

```
CMM Create MailMan Message          LR List Reminder Definitions
DFE Delete Exchange File Entry      LWH Load Web Host File
IFE Install Exchange File Entry     RI Reminder Definition Inquiry
IH Installation History              RP Repack
Select Action: Next Screen// lwh    Load Web Host File
Input the URL for the .prd file
https://[REDACTED]/UPDATE_2_0_370.PRD
```

At the **Select Action:** prompt, enter **LWH** for Load Web Host File

At the **Input the url for the .prd file:** prompt, type the following web address:

https://REDACTD/UPDATE 2 0 370.PRD

Follow the steps below to install

UPDATE_2_0_370 VA-TELERETINOPATHY UPDATES.

Search and locate an entry titled UPDATE_2_0_370 VA-TELERETINOPATHY UPDATE in reminder exchange.

```
512 UPDATE_2_0_370 VA-TELERETINOPA ←
    THY UPDATES
513 UPDATE_2_0_371 VA-WH BREAST
    CANCER RISK ASSESSMENT
514 UPDATE_2_0_372 VA-TELESTROKE
    DIALOG UPDATES
515 UPDATE_2_0_372 VA-TELESTROKE
    DIALOG UPDATES
516 UPDATE_2_0_39 VA-OSP SUICIDE
    SAFETY PLAN
+      + Next Screen  - Prev Screen  ?? More Actions  >>>
CFE Create Exchange File Entry      LHF Load Host File
CHF Create Host File                LMM Load MailMan Message
CMM Create MailMan Message          LR List Reminder Definitions
DFE Delete Exchange File Entry      LWH Load Web Host File
IFE Install Exchange File Entry     RI Reminder Definition Inquiry
IH Installation History              RP Repack
Select Action: Next Screen// ife    Install Exchange File Entry
Enter a list or range of numbers (1-751): 512
```

At the **Select Action** prompt, enter **IFE** for Install Exchange File Entry

Enter the number that corresponds with your entry **UPDATE_2_0_370 VA-TELERETINOPATHY UPDATE** (in this example it is entry 509 it will vary by site). The date of the exchange file should be 10/024/2023.

Description:

The following Clinical Reminder items were selected for packing:

REMINDER DIALOG

- VA-TELERETINOPATHY READER NOTE
- VA-TELERETINOPATHY IMAGER NOTE
- VA-EYE CARE AT-RISK SCREEN (NURSING)
- VA-EYE CARE AT-RISK SCREEN (PROV)
- VA-GP TECS TECH-DIABETES
- VA-TELE-EYE RESULTS LETTER

REMINDER DEFINITION

- VA-EYE CARE AT-RISK SCREEN (NURSE)

```
+      Enter ?? for more actions      >>>
IA  Install all Components              IS  Install Selected Component
Select Action: Next Screen// IA  Install all Components  ....
```

At the **Select Action** prompt, type **IA** for Install all Components and hit enter.

Select Action: Next Screen// **IA Install all Components**

You will see several prompts, for all new entries you will choose **I to Install**

- **Overwrite** existing health factors, reminder taxonomies, reminder definitions, and health summary types.
- **Skip** all exams that already exists but the packed component is different.
- **Merge** existing reminder terms.

Six dialogs need to be installed. Follow these steps to install each one.

You will be prompted to install the **first** reminder dialog component – **VA-EYE CARE AT-RISK SCREEN (PROV)**

```
Packed reminder dialog: VA-EYE CARE AT-RISK SCREEN (PROV)
```

Item	Seq.	Dialog Findings	Type	Exists
1		VA-EYE CARE AT-RISK SCREEN (PROV)	dialog	X
2	5	VA-EYE CARE AT-RISK SCREEN PROVIDER NOTE VERSION Finding: *NONE*	element	X
3	10	VA-GP EYE CARE AT-RISK SCREEN PROVIDER-REFUSED * Finding: *NONE*	group	X
4	10.5	VA-GP EYE CARE AT-RISK PATIENT HISTORY Finding: *NONE*	group	X
5	10.5.5	VA-OBJ EYE CARE AT-RISK PT HX-DIABETES Finding: *NONE*	element	X
6	10.5.10	VA-OBJ EYE CARE AT-RISK PT HX-MACULAR DEGENERATION Finding: *NONE*	element	X
7	10.5.15	VA-OBJ EYE CARE AT-RISK PT HX-GLAUCOMA	element	X

```
+      + Next Screen  - Prev Screen  ?? More Actions
DD  Dialog Details      DT  Dialog Text      IS  Install Selected
DF  Dialog Findings    DU  Dialog Usage    QU  Quit
DS  Dialog Summary     IA  Install All
Select Action: Next Screen// IA  Install All
```

At the **Select Action** prompt, type **IA** to install the dialog – **VA-EYE CARE AT-RISK SCREEN (PROV)**

Select Action: Next Screen// **IA Install All**

Install reminder dialog and all components with no further changes: Y// **Yes**

```
Packed reminder dialog: VA-EYE CARE AT-RISK SCREEN (PROV)
VA-EYE CARE AT-RISK SCREEN (PROV) (reminder dialog) installed from exchange file
```

Item	Seq.	Dialog Findings	Type	Exists
1		VA-EYE CARE AT-RISK SCREEN (PROV)	dialog	X
2	5	VA-EYE CARE AT-RISK SCREEN PROVIDER NOTE VERSION Finding: *NONE*	element	X
3	10	VA-GP EYE CARE AT-RISK SCREEN PROVIDER-REFUSED * Finding: *NONE*	group	X
4	10.5	VA-GP EYE CARE AT-RISK PATIENT HISTORY Finding: *NONE*	group	X
5	10.5.5	VA-OBJ EYE CARE AT-RISK PT HX-DIABETES Finding: *NONE*	element	X
6	10.5.10	VA-OBJ EYE CARE AT-RISK PT HX-MACULAR DEGENERATION Finding: *NONE*	element	X
7	10.5.15	VA-OBJ EYE CARE AT-RISK PT HX-GLAUCOMA	element	X

```
+      + Next Screen  - Prev Screen  ?? More Actions
DD  Dialog Details      DT  Dialog Text      IS  Install Selected
DF  Dialog Findings    DU  Dialog Usage    QU  Quit
DS  Dialog Summary     IA  Install All
Select Action: Next Screen// Q  Quit
```

After completing this dialog install, you will type Q.

You will then be prompted to install the **second** reminder dialog component –

VA-EYE CARE AT-RISK SCREEN (NURSING)

```
Packed reminder dialog: VA-EYE CARE AT-RISK SCREEN (NURSING)
```

Item	Seq.	Dialog Findings	Type	Exists
1		VA-EYE CARE AT-RISK SCREEN (NURSING)	dialog	X
2	2	VA-EYE CARE AT-RISK SCREEN NOTE VERSION Finding: *NONE*	element	X
3	5	VA-GP EYE CARE AT-RISK SCREEN Finding: *NONE*	group	X
4	5.5	VA-GP EYE CARE AT-RISK PATIENT HISTORY Finding: *NONE*	group	X
5	5.10	VA-BLANK LINE DIALOG & NOTE Finding: *NONE*	element	X
6	5.15	VA-OBJ TRI SCR N EXCLUSIONS Finding: *NONE*	element	X
7	5.20	VA-GP EYE CARE AT-RISK SCREENING-RESOLUTION Finding: *NONE*	group	X

```
+      + Next Screen  - Prev Screen  ?? More Actions
DD  Dialog Details      DT  Dialog Text      IS  Install Selected
DF  Dialog Findings    DU  Dialog Usage    QU  Quit
DS  Dialog Summary     IA  Install All
Select Action: Next Screen// IA  Install All
```

At the **Select Action** prompt, type **IA** to install the dialog – **VA-EYE CARE AT-RISK SCREEN (NURSING)**

Select Action: Next Screen// **IA Install All**

Install reminder dialog and all components with no further changes: Y// **Yes**

```
Packed reminder dialog: VA-EYE CARE AT-RISK SCREEN (NURSING)
VA-EYE CARE AT-RISK SCREEN (NURSING) (reminder dialog) installed from exchange f
```

Item	Seq.	Dialog Findings	Type	Exists
1		VA-EYE CARE AT-RISK SCREEN (NURSING)	dialog	X
2	2	VA-EYE CARE AT-RISK SCREEN NOTE VERSION Finding: *NONE*	element	X
3	5	VA-GP EYE CARE AT-RISK SCREEN Finding: *NONE*	group	X
4	5.5	VA-GP EYE CARE AT-RISK PATIENT HISTORY Finding: *NONE*	group	X
5	5.10	VA-BLANK LINE DIALOG & NOTE Finding: *NONE*	element	X
6	5.15	VA-OBJ TRI SCRNM EXCLUSIONS Finding: *NONE*	element	X
7	5.20	VA-GP EYE CARE AT-RISK SCREENING-RESOLUTION Finding: *NONE*	group	X

```
+      + Next Screen  - Prev Screen  ?? More Actions
DD  Dialog Details      DT  Dialog Text      IS  Install Selected
DF  Dialog Findings    DU  Dialog Usage    QU  Quit
DS  Dialog Summary     IA  Install All
Select Action: Next Screen// Q  Quit
```

After completing this dialog install, you will type Q.

You will then be prompted to install the **third** reminder dialog component – **VA-TELE-EYE RESULTS LETTER**

```
Packed reminder dialog: VA-TELE-EYE RESULTS LETTER [NATIONAL DIALOG]
```

Item	Seq.	Dialog Findings	Type	Exists
1		VA-TELE-EYE RESULTS LETTER	dialog	
2	5	VA-TXT TELE-EYE RESULT LETTER VERSION Finding: *NONE*	element	
3	10	VAL-GP TELE-EYE RESULTS LETTER Finding: VA-TRR RESULTS NOTIFICATION COMPLETED (HEALTH FACTOR)	group	X
4	10.5	VA-TXT TELE-EYE RESULT LETTER-RESULTS/RECOMMEND ATIONS Finding: *NONE*	element	
5	10.10	VAL-TXT TELE-EYE RESULT LETTER-LOCAL INFO Finding: *NONE*	element	
6	10.15	VA-TXT TELE-EYE RESULT LETTER-REPORT/TOBACCO Finding: *NONE*	element	

```
+      + Next Screen  - Prev Screen  ?? More Actions
DD  Dialog Details      DT  Dialog Text      IS  Install Selected
DF  Dialog Findings    DU  Dialog Usage    QU  Quit
DS  Dialog Summary     IA  Install All
Select Action: Next Screen// IA  Install All
```

At the **Select Action** prompt, type **IA** to install the dialog – **VA-TELE-EYE RESULTS LETTER**

Select Action: Next Screen// **IA Install All**

Install reminder dialog and all components with no further changes: Y// **Yes**

```

Packed reminder dialog: VA-TELE-EYE RESULTS LETTER [NATIONAL DIALOG]
VA-TELE-EYE RESULTS LETTER (reminder dialog) installed from exchange file.

```

Item	Seq.	Dialog Findings	Type	Exists
1		VA-TELE-EYE RESULTS LETTER	dialog	X
2	5	VA-TXT TELE-EYE RESULT LETTER VERSION Finding: *NONE*	element	X
3	10	VAL-GP TELE-EYE RESULTS LETTER Finding: VA-TRR RESULTS NOTIFICATION COMPLETED (HEALTH FACTOR)	group	X X
4	10.5	VA-TXT TELE-EYE RESULT LETTER-RESULTS/RECOMMEND ATIONS Finding: *NONE*	element	X
5	10.10	VAL-TXT TELE-EYE RESULT LETTER-LOCAL INFO Finding: *NONE*	element	X
6	10.15	VA-TXT TELE-EYE RESULT LETTER-REPORT/TOBACCO Finding: *NONE*	element	X

```

+      + Next Screen  - Prev Screen  ?? More Actions
DD  Dialog Details      DT  Dialog Text      IS  Install Selected
DF  Dialog Findings    DU  Dialog Usage    QU  Quit
DS  Dialog Summary     IA  Install All
Select Action: Next Screen// Q  Quit

```

After completing this dialog install, you will type Q.

You will then be prompted to install the **fourth** reminder dialog component - **VA-GP TECS TECH-DIABETES**

```

Packed reminder dialog: VA-GP TECS TECH-DIABETES [NATIONAL DIALOG]

```

Item	Seq.	Dialog Findings	Type	Exists
1		VA-GP TECS TECH-DIABETES* Finding: *NONE*	group	X
2	5	VA-TXT TECS TECH DIABETES-NO Finding: *NONE*	element	X
3	10	VA-GP TECS TECH DIABETES-YES Finding: VA-TECS EYE DX HX-DIABETES (HEALTH FACTOR)	group	X X
4	10.5	VA-TXT TECS TECH DIABETES-DURATION Finding: *NONE*	element	X
5	10.10	VA-GP TECS TECH DIABETES-INSULIN Finding: *NONE*	group	X
6	10.10.5	VA-TXT TECS TECH INSULIN-NO Finding: *NONE*	element	X
7	10.10.10	VA-HF TECS TECH INSULIN-YES	element	X

```

+      + Next Screen  - Prev Screen  ?? More Actions
DD  Dialog Details      DT  Dialog Text      IS  Install Selected
DF  Dialog Findings    DU  Dialog Usage    QU  Quit
DS  Dialog Summary     IA  Install All
Select Action: Next Screen// IA  Install All

```

At the **Select Action** prompt, type **IA** to install the dialog - **VA-GP TECS TECH-DIABETES**

Select Action: Next Screen// **IA Install All**

Install reminder dialog and all components with no further changes: Y// **Yes**

Packed reminder dialog: VA-GP TECS TECH-DIABETES [NATIONAL DIALOG]
 VA-GP TECS TECH-DIABETES (group) installed from exchange file.

Item	Seq.	Dialog Findings	Type	Exists		
1		VA-GP TECS TECH-DIABETES* Finding: *NONE*	group	X		
2	5	VA-TXT TECS TECH DIABETES-NO Finding: *NONE*	element	X		
3	10	VA-GP TECS TECH DIABETES-YES Finding: VA-TECS EYE DX HX-DIABETES (HEALTH FACTOR)	group	X		
4	10.5	VA-TXT TECS TECH DIABETES-DURATION Finding: *NONE*	element	X		
5	10.10	VA-GP TECS TECH DIABETES-INSULIN Finding: *NONE*	group	X		
6	10.10.5	VA-TXT TECS TECH INSULIN-NO Finding: *NONE*	element	X		
7	10.10.10	VA-HF TECS TECH INSULIN-YES	element	X		
+ + Next Screen - Prev Screen ?? More Actions						
DD	Dialog Details		DT	Dialog Text	IS	Install Selected
DF	Dialog Findings		DU	Dialog Usage	QU	Quit
DS	Dialog Summary		IA	Install All		
Select Action: Next Screen// Q Quit						

After completing this dialog install, you will type **Q**.

You will then be prompted to install the **fifth** reminder dialog – VA-TELERETINOPATHY READER NOTE

Packed reminder dialog: VA-TELERETINOPATHY READER NOTE [NATIONAL DIALOG]

Item	Seq.	Dialog Findings	Type	Exists		
1		VA-TELERETINOPATHY READER NOTE	dialog	X		
2	5	VA-TELE-EYE READER NOTE VERSION Finding: *NONE*	element	X		
3	10	VA-GP TRR TELE-EYE CARE IMAGE QUALITY Finding: *NONE*	group	X		
4	10.5	VA-HF TRR TELE-EYE EXAM QUAL ADEQUATE Finding: VA-TRR IMAGE ADEQUATE ICD CODE (REMINDER TAXONOMY) Add. Finding: VA-TRR IMAGE QUAL ADEQUATE (HEALTH FACTOR)	element	X		
5		FORCED PRIMARY DIAGNOSIS	forced	X		
6		PXRM COMMENT	prompt	X		
7	10.10	VA-GP TRR TELE-EYE EXAM QUAL INADEQUATE Finding: VA-TRR IMAGE INADEQUATE ICD CODE (REMINDER TAXONOMY)	group	X		
+ + Next Screen - Prev Screen ?? More Actions						
DD	Dialog Details		DT	Dialog Text	IS	Install Selected
DF	Dialog Findings		DU	Dialog Usage	QU	Quit
DS	Dialog Summary		IA	Install All		
Select Action: Next Screen// IA Install All						

At the **Select Action** prompt, type **IA** to install the dialog – VA-TELERETINOPATHY READER NOTE

Select Action: Next Screen// **IA** Install All

Install reminder dialog and all components with no further changes: Y// **Yes**

```

Packed reminder dialog: VA-TELERETINOPATHY READER NOTE [NATIONAL DIALOG]
VA-TELERETINOPATHY READER NOTE (reminder dialog) installed from exchange file.
Item  Seq. Dialog Findings                                     Type  Exists
  1    VA-TELERETINOPATHY READER NOTE                         dialog  X
  2    5 VA-TELE-EYE READER NOTE VERSION                       element  X
      Finding: *NONE*
  3    10 VA-GP TRR TELE-EYE CARE IMAGE QUALITY                group   X
      Finding: *NONE*
  4    10.5 VA-HF TRR TELE-EYE EXAM QUAL ADEQUATE              element  X
      Finding: VA-TRR IMAGE ADEQUATE ICD CODE (REMINDER TAXONOMY)  X
      Add. Finding: VA-TRR IMAGE QUAL ADEQUATE (HEALTH FACTOR)    X
  5    FORCED PRIMARY DIAGNOSIS                                forced   X
  6    PXR COMMENT                                             prompt   X
  7    10.10 VA-GP TRR TELE-EYE EXAM QUAL INADEQUATE          group   X
      Finding: VA-TRR IMAGE INADEQUATE ICD CODE (REMINDER
      TAXONOMY)                                                X
+      + Next Screen  - Prev Screen  ?? More Actions
DD  Dialog Details      DT  Dialog Text      IS  Install Selected
DF  Dialog Findings     DU  Dialog Usage     QU  Quit
DS  Dialog Summary      IA  Install All
Select Action: Next Screen// Q  Quit

```

After completing this dialog install, you will type **Q**.

You will then be prompted to install the **sixth** reminder dialog – **VA-TELERETINOPATHY IMAGER NOTE**

```

Packed reminder dialog: VA-TELERETINOPATHY IMAGER NOTE [NATIONAL DIALOG]
Item  Seq. Dialog Findings                                     Type  Exists
  1    VA-TELERETINOPATHY IMAGER NOTE                         dialog  X
  2    5 VA-TELE-EYE IMAGER NOTE VERSION                       element  X
      Finding: *NONE*
  3    10 VA-GP TRI IMAGER CONSULT                             group   X
      Finding: *NONE*
  4    10.5 VA-GP TRI IMAGING CONSULT INITIAL                  group   X
      Finding: *NONE*
  5    10.5.5 VA-TXT TRI IMAGING DEMOGRAPHICS                  element  X
      Finding: *NONE*
  6    10.5.10 VA-GP TRI IMAGER'S CONSULT INITIAL RESULT      group   X
      Finding: *NONE*
  7    10.5.10.5 VA-TXT TRI BLANK SPACE                        element  X
      Finding: *NONE*
+      + Next Screen  - Prev Screen  ?? More Actions
DD  Dialog Details      DT  Dialog Text      IS  Install Selected
DF  Dialog Findings     DU  Dialog Usage     QU  Quit
DS  Dialog Summary      IA  Install All
Select Action: Next Screen// IA  Install All

```

At the **Select Action** prompt, type **IA** to install the dialog – **VA-TELERETINOPATHY IMAGER NOTE**

Select Action: Next Screen// **IA Install All**

Install reminder dialog and all components with no further changes: Y// **Yes**

```

Packed reminder dialog: VA-TELERETINOPATHY IMAGER NOTE [NATIONAL DIALOG]
VA-TELERETINOPATHY IMAGER NOTE (reminder dialog) installed from exchange file.
Item  Seq. Dialog Findings                                     Type  Exists
  1    VA-TELERETINOPATHY IMAGER NOTE                         dialog  X
  2    5 VA-TELE-EYE IMAGER NOTE VERSION                       element  X
      Finding: *NONE*
  3    10 VA-GP TRI IMAGER CONSULT                             group   X
      Finding: *NONE*
  4    10.5 VA-GP TRI IMAGING CONSULT INITIAL                  group   X
      Finding: *NONE*
  5    10.5.5 VA-TXT TRI IMAGING DEMOGRAPHICS                 element  X
      Finding: *NONE*
  6    10.5.10 VA-GP TRI IMAGER'S CONSULT INITIAL RESULT      group   X
      Finding: *NONE*
  7    10.5.10.5 VA-TXT TRI BLANK SPACE                       element  X
      Finding: *NONE*
+      + Next Screen  - Prev Screen  ?? More Actions
DF  Dialog Findings  DU  Dialog Usage  QU  Quit
DS  Dialog Summary  IA  Install All
Select Action: Next Screen// Q  Quit

```

After completing this dialog install, you will type **Q**.

```

Description:
The following Clinical Reminder items were selected for packing:
REMINDER DIALOG
  VA-TELERETINOPATHY READER NOTE
  VA-TELERETINOPATHY IMAGER NOTE
  VA-EYE CARE AT-RISK SCREEN (NURSING)
  VA-EYE CARE AT-RISK SCREEN (PROV)
  VA-GP TECS TECH-DIABETES
  VA-TELE-EYE RESULTS LETTER

REMINDER DEFINITION
  VA-EYE CARE AT-RISK SCREEN (NURSE)
+      + Next Screen  - Prev Screen  ?? More Actions  >>>
IA  Install all Components  IS  Install Selected Component
Select Action: Next Screen// Q  Q

```

You will then be returned to this screen. At the Select Action prompt, type **Q**.

Install complete.

Post Installation

1. Confirm all local findings are still mapped to the following reminder terms:

- a. TELE EYE GLAUCOMA EXAM-ABNORMAL
- b. TELE EYE GLAUCOMA EXAM-NORMAL

- c. TELE EYE GLAUCOMA EXAM-UNABLE TO ASSESS
- d. TELE EYE MACULA EXAM-ABNORMAL
- e. TELE EYE MACULA EXAM-NORMAL
- f. TELE EYE MACULA EXAM-UNABLE TO ASSESS
- g. TELE EYE RETINOPATHY EXAM-ABNORMAL
- h. TELE EYE RETINOPATHY EXAM-NORMAL
- i. TELE EYE RETINOPATHY EXAM-UNABLE TO ASSESS
- j. TELE EYE SCREEN ORDERS
- k. VA-OUTSIDE EYE AT-RISK SCREEN
- l. VA-EYE AT RISK NSG SCREEN REFUSAL
- m. VA-EYE AT RISK NSG SCREEN INCORRECT DM DIAGNOSIS
- n. VA-BL TELE EYE INCORRECT DIABETIC DX

2. Add orders as an additional finding item to the following reminder dialog elements.

- a. Element VAL-OI EYE CARE AT-RISK RETINAL IMAGE CONSULT
 - i. Add local quick order for Retinal Image Consult

NAME: **VAL-OI EYE CARE AT-RISK RETINAL IMAGE CONSULT**

DISABLE:

CLASS: LOCAL//

SPONSOR:

REVIEW DATE:

RESOLUTION TYPE: ORDERED//

ORDERABLE ITEM:

FINDING ITEM: Enter "Q.", then name of local quick order for Retinal Image Consult

- b. Element VAL-OI EYE CARE AT-RISK EYE GLASS CONSULT
 - i. Add local quick order/process for scheduling an eye clinic appointment for glasses

NAME: **VAL-OI EYE CARE AT-RISK EYE GLASS CONSULT**

DISABLE:

CLASS: LOCAL//

SPONSOR:

REVIEW DATE:

RESOLUTION TYPE: ORDERED//

ORDERABLE ITEM:

FINDING ITEM: Enter "Q.", then name of local quick order for scheduling an eye clinic appointment for glasses

- c. Element VAL-OI EYE CARE AT-RISK COMP EYE EXAM CONSULT
 - i. Add local quick order/process for scheduling an eye clinic appointment

NAME: **VAL-OI EYE CARE AT-RISK COMP EYE EXAM CONSULT**

DISABLE:

CLASS: LOCAL//

SPONSOR:

REVIEW DATE:

RESOLUTION TYPE: ORDERED//

ORDERABLE ITEM:

FINDING ITEM: Enter "Q.", then name of local quick order for scheduling an eye clinic appointment

NAME: **VA-TXT TRR RECOMMENDATION OTH OPTIONS-NUTRITION REFERRAL**

FINDING ITEM: Enter "Q.", then name of local quick order for referral to Nutrition.

3. Map local hemoglobin A1C lab tests to the Health Summary Type TRI HGBA1C.

Begin in the Health Summary Coordinator's Menu

- 1 Print Health Summary Menu ...
- 2 Build Health Summary Type Menu ...
- 3 Set-up Batch Print Locations
- 4 List Batch Health Summary Locations
- 5 CPRS Reports Tab 'Health Summary Types List' Menu ...

Select Health Summary Coordinator's Menu <TEST ACCOUNT> Option: **2** Build Health Summary Type Menu

- 1 Create/Modify Health Summary Type
- 2 Delete Health Summary Type
- 3 Health Summary Objects Menu ...
- 4 Information Menu ...
- 5 Print Health Summary Menu ...

Select Build Health Summary Type Menu <TEST ACCOUNT> Option: **1** Create/Modify Health Summary Type

Select Health Summary Type: **Tri Hgba1c**

WARNING: You are about to edit a Health Summary Type that is being used by a Health Summary Object. Changing the structure of this Health Summary Type will alter how the Object will display.

Do want to continue? NO// **YES**

NAME: TRI HGBA1C//

TITLE: TRI HGBA1C//

SUPPRESS PRINT OF COMPONENTS WITHOUT DATA: yes//

SUPPRESS SENSITIVE PRINT DATA: NO SSN//

LOCK:

OWNER:

Do you wish to review the Summary Type structure before continuing? NO// **YES**

HEALTH SUMMARY TYPE INQUIRY

Type Name: TRI HGBA1C

Title: TRI HGBA1C

Owner:

SUPPRESS PRINT OF COMPONENTS WITHOUT DATA: yes

SUPPRESS SENSITIVE PRINT DATA: NO SSN

Abb	Ord	Component Name	Max Occ	Time	Hos Loc	ICD Text	Pro Nar	CPT Mod	Selection
SLT	5	Lab Tests Selected	3	2Y	Confirm the max occurrence and time limit matches the highlighted text				

* = Disabled Components

Select COMPONENT: **SLT** LAB TESTS SELECTED SLT

LAB TESTS SELECTED is already a component of this summary.

Select one of the following:

- E Edit component parameters
- D Delete component from summary

Select Action: **Edit component parameters**

SUMMARY ORDER: 5// 5

OCCURRENCE LIMIT: **3**

TIME LIMIT: **2Y**

HEADER NAME: Lab Tests Selected//

No selection items chosen.

Select new items one at a time in the sequence you want them displayed.

You may select up to 99 items.

Select SELECTION ITEM: **Enter local lab(s) for hemoglobin A1C**

4. Map local drug(s) for Semaglutide to the VA-BL SEMAGLUTIDE RX reminder term.

Select Reminder Term: **VA-BL SEMAGLUTIDE RX**

Choose from:

FINDING ITEM: **SEMAGLUTIDE** **Enter all local drug(s) for Semaglutide**

BEGINNING DATE/TIME: **T-18M**

ENDING DATE/TIME:

OCCURRENCE COUNT:

RXTYPE: **A**

USE START DATE: **YES**

CONDITION:
CONDITION CASE SENSITIVE: NO//
USE STATUS/COND IN SEARCH: **YES**

Statuses already defined for this finding item: **ACTIVE**

Select one of the following:

- A ADD STATUS
- D DELETE A STATUS
- S SAVE AND QUIT
- Q QUIT WITHOUT SAVING CHANGES

5. Add local data objects to the VAL-GP TELE-EYE RESULTS LETTER group.

NAME: **VAL-GP TELE-EYE RESULTS LETTER**
DISABLE:
CLASS: LOCAL//
SPONSOR: OFFICE OF TELEHEALTH SERVICES STORE AND FORWARD TELEHEALTH//
REVIEW DATE:
RESOLUTION TYPE:
ORDERABLE ITEM:
Finding item: HF VA-TRR RESULTS NOTIFICATION COMPLETED
FINDING ITEM: VA-TRR RESULTS NOTIFICATION COMPLETED//

Additional findings: none
Select ADDITIONAL FINDING:
GROUP CAPTION:
PUT A BOX AROUND THE GROUP:
SHARE COMMON PROMPTS:
MULTIPLE SELECTION:
HIDE/SHOW GROUP: SHOW//
GROUP HEADER DIALOG TEXT:. . .

. . .
Placeholder for PATIENT NAME\<\
Placeholder for PATIENT STREET\<\
Placeholder for PATIENT CITY\<\
Placeholder for PATIENT STATE\<\
Placeholder for PATIENT ZIP CODE

Dear Veteran Placeholder for Patient's First Name Patient's Last Name:

Edit? NO// **YES**

==[WRAP]==[INSERT]===== < GROUP HEADER DIALOG TEX[Press <PF1>H for help]=====

Placeholder for Today's Date

**Placeholder for PATIENT NAME\<\
Placeholder for PATIENT STREET\<\
Placeholder for PATIENT CITY\<**

Placeholder for PATIENT STATE\
Placeholder for PATIENT ZIP CODE

Dear Veteran Placeholder for Patient's First Name Patient's Last Name:

Replace highlighted text with local objects.

<=====T=====T=====T=====T=====T=====T=====T=====T=====T>=====

6. Add local data object to the VAL-OBJ MOST RECENT EYE EXAM DATES element. The local object should display the most recent eye exam dates. Remove placeholder text if no local objects exists.

NAME: VAL-OBJ MOST RECENT EYE EXAM DATES

DISABLE:

CLASS: LOCAL//

SPONSOR: OFFICE OF TELEHEALTH SERVICES STORE AND FORWARD TELEHEALTH//

REVIEW DATE:

RESOLUTION TYPE:

ORDERABLE ITEM:

Finding item:

Additional findings: none

Select ADDITIONAL FINDING:

SHARE COMMON PROMPTS:

GROUP HEADER DIALOG TEXT:. . .

Last Retinal Evaluation:

Place local data object for most recent eye exam dates

Edit? NO// YES

==[WRAP]==[INSERT]===== < GROUP HEADER DIALOG TEX[Press <PF1>H for help]=====

Last Retinal Evaluation:

Place local data object for most recent eye exam dates

Replace highlighted text with local objects or remove it from group if no object exists.

<=====T=====T=====T=====T=====T=====T=====T=====T=====T>=====

7. Element VAL-TXT TELE-EYE RESULT LETTER-LOCAL INFO is a local class in the VA-TELE-EYE RESULTS LETTER dialog, so sites can revise text to meet local processes.

NAME: VAL-TXT TELE-EYE RESULT LETTER-LOCAL INFO

DISABLE:

CLASS: LOCAL//

SPONSOR:

REVIEW DATE:

RESOLUTION TYPE: OTHER//

ORDERABLE ITEM:

FINDING ITEM:

Additional findings: none

Select ADDITIONAL FINDING:

DIALOG/PROGRESS NOTE TEXT:

If you have not received a call within 30 days of the recommended appointment date, or have vision or eye problems prior to your next appointment, please call to schedule an appointment:

Scheduling Number: XXX-XXX-XXXX

Edit? NO// **YES**

==[WRAP]==[INSERT]===== < DIALOG/PROGRESS NOTE TE[Press <PF1>H for help]=====

If you have not received a call within 30 days of the recommended appointment date, or have vision or eye problems prior to your next appointment, please call to schedule an appointment:

Scheduling Number: XXX-XXX-XXXX

Replace highlighted text as needed. Add scheduling phone number.

8. Make the Dialog template available to be attached to a Personal or Shared template or to be able to attach the dialog to a progress note title.

Select Reminder Managers Menu <TEST ACCOUNT> Option: **CP**

- CA Add/Edit Reminder Categories
- CL CPRS Lookup Categories
- CS CPRS Cover Sheet Reminder List
- MH Mental Health Dialogs Active
- PN Progress Note Headers
- RA Reminder GUI Resolution Active
- TIU TIU Template Reminder Dialog Parameter**
- DEVL Evaluate Coversheet List on Dialog Finish
- DL Default Outside Location
- PT Position Reminder Text at Cursor
- LINK Link Reminder Dialog to Template
- TEST CPRS Coversheet Time Test
- NP New Reminder Parameters
- GEC GEC Status Check Active
- WH WH Print Now Active

Select CPRS Reminder Configuration <TEST ACCOUNT> Option: **TIU TIU Template Reminder Dialog Parameter**

Reminder Dialogs allowed as Templates may be set for the following:

- 1 User USR [choose from NEW PERSON]
- 3 Service SRV [choose from SERVICE/SECTION]
- 4 Division DIV [SALT LAKE CITY]

5 System SYS [NATREM.FO-SLC.MED.VA.GOV]

Enter selection: **5** Choose the appropriate number for your site. Your site may do this by System or other levels. The example below uses SYSTEM level.

Setting Reminder Dialogs allowed as Templates for System: NATREM.FO-SLC.MED.VA.GOV
Select Display Sequence: **?**

388 VA-COVID-19 TELEPHONE FOLLOW-UP
389 VA-SRA-CONSULT RESPONSE (D)

When you type a question mark above, you will see the list of #'s (with dialogs) that are already taken. Choose a number **NOT** on this list. For this example, 390 and 391 are not present so those numbers can be used.

Select Display Sequence: **390**
Are you adding 390 as a new Display Sequence? Yes// **YES**

Display Sequence: 390// **<Enter>**
Clinical Reminder Dialog: **VA-TELE-EYE RESULTS LETTER**
reminder dialog NATIONAL
...OK? Yes// **Yes**

9. Setup of Note Title.

The note title to create to associate with this template is: **TELE-EYE RESULTS LETTER**


VHA Enterprise Standard Title: **LETTERS**

You can add this to the appropriate document class.

Note: Comprehensive information on Creating new Document Classes and Titles can be found beginning on page 54 of the TIU/ASU Implementation Guide at the following link:
[http://www.va.gov/vdl/documents/Clinical/CPRS-Text_Integration_UTILITY_\(TIU\)/tiuim.pdf](http://www.va.gov/vdl/documents/Clinical/CPRS-Text_Integration_UTILITY_(TIU)/tiuim.pdf)

10. Associate the reminder dialog with the note title in CPRS

In CPRS, open the Template Editor, expand Document Titles, select the appropriate folder (if applicable), then click New Template (top right-hand corner).

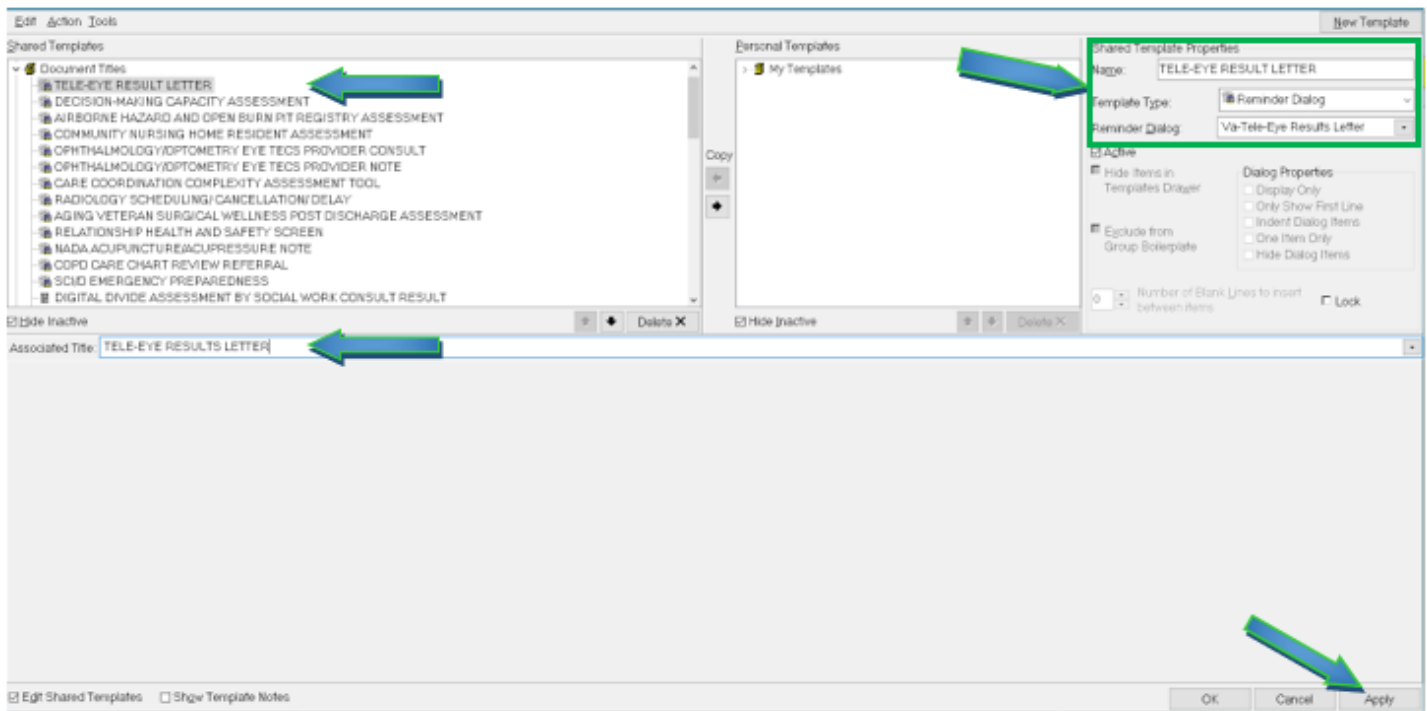

Type in the following information:

Name: **TELE-EYE RESULTS LETTER**

Template Type: **Reminder Dialog**

Reminder Dialog: **VA-TELE-EYE RESULTS LETTER**

Next, type in the Associated Title name -- then click Apply.